

Lifelong

## Lifelong Mini Gun Massager



4 Massage  
Heads



No wire  
required



3200  
Strokes/min

**1 YEAR  
WARRANTY**

## INTRODUCTION:

Welcome to the Lifelong family,

Dear Customer,

Congratulations on purchasing Lifelong Mini Gun Massager and a warm welcome to the Lifelong family.

At Lifelong, we take pride in developing quality products for your home, health and happiness. Through constant innovation in our products, we have been able to develop a wide spectrum of products across Home, Kitchen, Grooming and Wellness categories. We strive to enrich the lives of the modern consumer every day. Keeping that in mind, we have offered all our products at a fair and honest price.

We hope that the Mini Gun Massager serves you and your family's needs without any hassles.

Thanks again for becoming a part of the Lifelong family. Enjoy using the product!

Warm Regards,

Team Lifelong

## PRODUCT USE:


### On / off

- Press and hold the start key for three seconds to start.
- Press and hold the start key for three seconds to shut down.


### Force adjustment

- The starting gear is first gear.
- Click the start key once to increase one gear.
- The maximum gear is fourth gear.
- Click the start key once in fourth gear to change to first gear.


Note: in order to protect bones and muscles, avoid long-term fixed-point massage, resulting in damage, after 10 minutes it automatically stops with the intelligent protection mode.

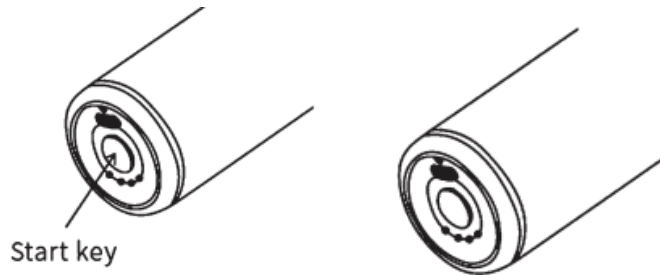

## PRODUCT OVERVIEW:

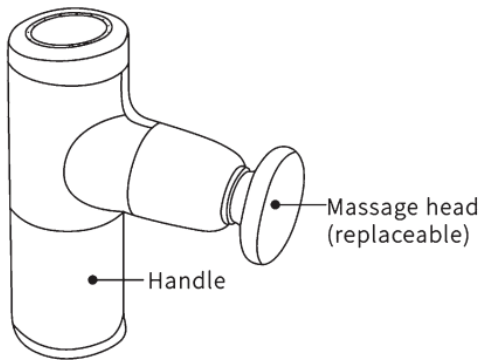


## MASSAGE HEAD:

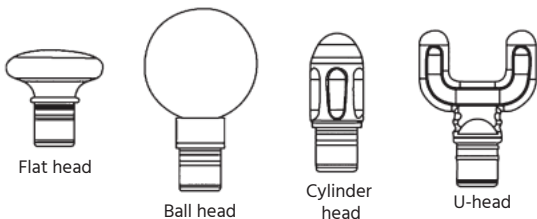

- Flat Head: It is suitable for relaxation and shaping of various muscle parts.
- Ball Head: Suitable for massaging all muscle parts, such as palms, soles etc.
- Cylinder Head: It is applicable to impact deep tissues, such as meridian joints etc.
- U-Head: Apply to massage the neck and achilles tendon of spine.

## CHANGE MESSAGE HEAD:

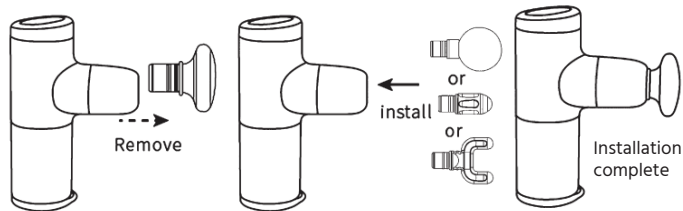


## INDICATOR LIGHT:


- The indicator light indicates the gear.
- The first gear is when one light is on, and the second gear is when two lights are on.
- Third gear when three lights are on, fourth gear when four lights are on.

## SIGNAL LIGHT:

- The signal light is always green when the machine is turned on.
- Four indicator lights are on once in turn.
- The signal light is red when charging.
- The signal light is green when the battery is fully charged.
- When the power is turned on, the signal lamp with sufficient power flashes 3 times every 5 seconds, 2 times every 5 seconds when the power is lower than 60%, and 1 time every 5 seconds when the power is lower than 30%.

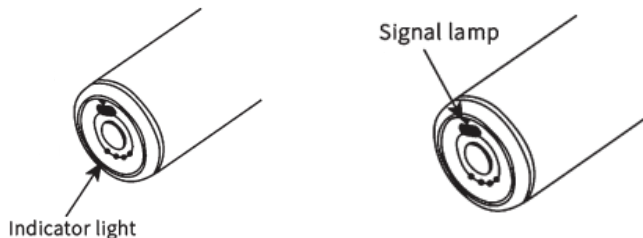

## CHARGE:

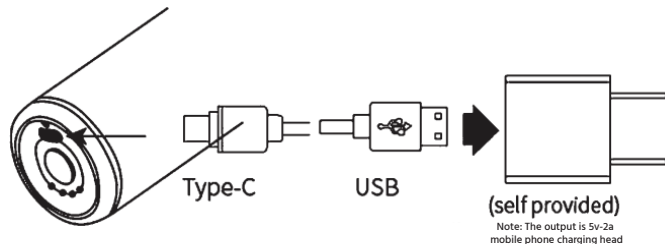

- Please be sure to fully charge before the first use.
- The signal light is red when charging, and green after charging.
- The charging time is about 4h, and the maximum usage range is about 11.5h.
- When the power of the mini gun is less than 10%, the signal light flashes red.
- When it turns red and charge as soon as possible.

**\*To ensure safety, do not use the mini gun when charging.**

**\*In order to protect the battery life, the signal light flashes red and indicates to charge it as soon as possible.**

## BENEFITS:

- Pain and spasm due to muscle contusion, sprain and strain.
- To help the swelling flow.
- Loose and thickened connective tissue and fascia.
- Reduce the accumulation of lactic acid in muscle.
- Increase joint mobility.
- Eliminate muscle fatigue.

## DO NOT USE:

- Aneurysm, bleeding, blood thinners.
- Heart disease, with pacemaker or defibrillator.
- Pregnancy, cancer.
- The device cannot be used within 3 inches of internal fixation within 90 days after internal fixation.
- Sensitive area head, face, cervical vertebra, spine and vertebral body.
- The position of tibia, dorsum of foot, dorsum of hand, etc. near the skeleton with less muscle coverage.

## PRODUCT PARAMETERS:

- Rated voltage: 8.4v
- Rated current: 2A
- Input parameter: 5V-2A
- Maximum speed: 3200 strokes/ min
- Battery capacity: 2000mAh

## Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log on to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.

### Terms and Conditions:

- This warranty is void if:
  - The completely filled warranty card is not presented at the time of servicing the product.
  - The product is not operated according to the instructions given under the user manual.
- Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
- Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
- Product label specifying the model number, serial number and production code has been removed and altered.
- Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
- Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
- Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
- Repairs or replacements will be carried out by authorized service provider.
- During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised service provider.
- All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
- Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
- Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
- However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
- No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

For queries and complaints:  
please contact: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)