

2min Timer The toothbrush will stop automatically after 2min brushing which is recommended by the dentists. How to charge:

Please move to the next section of your mouth every 30s

30s Interval Reminder

interval.

Before using for the first time, please charge the device. Connect the device to the USB charging cable. Insert the USB cable into any USB power adapter and connect it to a power supply.

The 3 lights will start flashing intermittently indicating that the brush is now charging. It takes about 3 hours to complete the charging. All the 3 lights will glow steadily once the device is completely charged. While charging, the device cannot

When battery power is lower than 20%, the 3 lights will flash at the same time. This serves as a reminder that the brush requires charging. During this period, the toothbrush will continue to function normally. When the battery power is almost depleted, the 3 lights will flash 5 times after which the device will shut down.

How to Use:

1) Attach the brush head and make sure it is locked correctly. 2) Wet the bristles and apply appropriate amount of toothpaste 3) Position the bristles at a 45 degree angle towards the



4) Turn the brush on when it is in your mouth to avoid splashing. 5) Gently guide and glide the brush along the gum-line

using only light pressure from your fingertips. Do NOT PRESS or SCRUB HARD. Slowly move the brush from tooth to tooth on the outer






6) Now clean the biting surfaces (top of the back teeth) by slowly moving across and letting the toothbrush do its job. FOR BEST RESULTS, MAKE SURE YOU CHANGE YOUR BRUSH HEAD EVERY 3 MONTHS. THIS CAN SIGNIFICANTLY REDUCE PLAQUE BUILDUP.

Changing the toothbrush head:

Changing the toothbrush head INSTALL: Slide the toothbrush head down on the handle. Turn clockwise until the bristles are aligned with the button.

DETACH: Turn the toothbrush head counterclockwise and gently pull out the brush head along the handle.

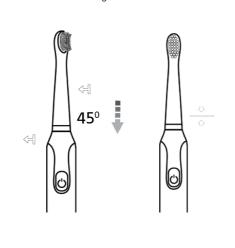

9

Consumer Warranty Card

Dear Customer.

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

All Lifelong consumer products are covered against manufacturing defects from the date of purchase.

Name of the product: _ Warranty coverage: 1 year Please note: Purchase reciept is necessary for warranty Customer Care: customercare@lifelongindia.com

Purchase Date:

Invoice Number:

Online Site: __

Address: _ Home Number Office Number: E-mail Address:

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days. 10

Terms and Conditions:

1. This warranty is void if: a. The completely filled warranty card is not presented at the

time of servicing the product. b. The product is not operated according to the instructions given under the user manual.

C. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.

d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service

organizations or persons. e. Product label specifying the model number, serial number

and production code has been removed and altered. f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.

g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received. h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-

household/non-domestic use. 2. Repairs or replacements will be carried out by authorized service provider.

3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its

4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except a set forth above, there are no other express or implied warranty and all warrantied, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted

by the lass. 6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.

7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law. 8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

12

Lifelong Online Retail Pvt Ltd Sth Floor, Unit No. 508, DLF Center, Saket, New Delhi -110017

 ${\it LifelongOnlineIn}~\ref{fig:poisson}$ **\$\square\$** +91-9711558877 Lifelong_Online 🖤 customercare@lifelongindia.com Lifelong.Online www.lifelongindiaonline.com

11