

Manufactured Market by:
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com

lifelong

Multi Effect Portable Smart Oral Irrigator



0.7mm Ultrafine
water flow



Detachable and
easy to clean



360 rotatable
nozzle



IPX7
waterproof

**1 YEAR
WARRANTY**

INTRODUCTION:

Dear Customer,

Congratulations!

We welcome you to the Lifelong family, as you welcome our Multi Effect Portable Smart Oral Irrigator to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Multi Effect Portable Smart Oral Irrigator. We know we can't!

Looking forward to building our bond for life, cheers!

Thank you,
Team Lifelong

TABLE OF CONTENTS:

About the product

Parts/accessories of product/ Marking of product

How to install

How to use

Device Cleaning


Care instructions

Safety for usage

Technical data/table for specifications

ABOUT THE PRODUCT:

The Lifelong Multi Effect Portable Smart Oral Irrigator gives you the best of a smooth water jet cleaning and pressure, effectively enhancing your oral health. It reaches places between the teeth and the gums that are hard to cover with traditional flossing. With a powerful 2000mAH battery and the 350 ml water tank, dental care has never been easier. The 360°rotatable nozzle gives you an all-rounded cleaning, allowing you to easily reach areas that are difficult to clean, and prevents tooth decay, dental plaque, dental calculus, gingival bleeding and dental hypersensitivity. Smile wider with Lifelong Portable Smart Oral Regulator. With more than 100 innovative products, no wonder millions of customers continue to choose Lifelong! It has been strongly trusted by loyal users for over five years. Made specially to keep you smiling brighter.

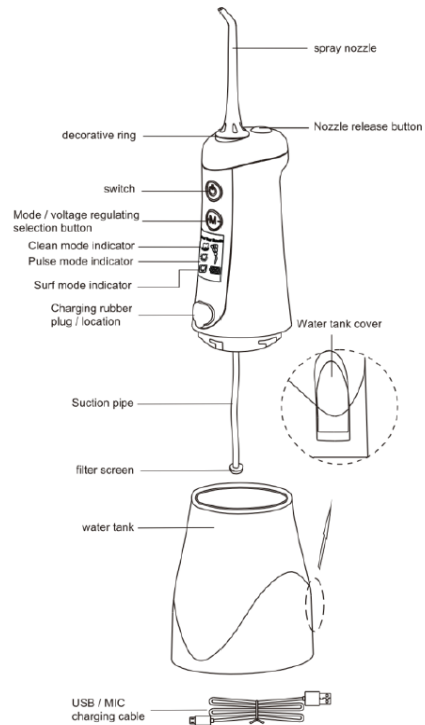

PRODUCT ACCESSORIES:

- 1 Lifelong Home Cordless Flosser
- 4 Classic Jet Tips
- 1 Jet Tip
- 1 Orthodontal Tip
- 1 Tongue Cleaner Tip
- 1 Periodontal Tip
- 1 USB cord for charging
- 1 User Manual

HOW TO INSTALL:

- 1 Insert the nozzle into the small opening on the top of the device. Always turn the flosser off before trying to replace the nozzle.
- 2 Gently press the nozzle release button to release the nozzle.
- 3 You can fill the water tank using two methods.
 - a Open the water tank cover, and hold the flosser horizontally before carefully filling water into the tank.
 - b Separate the water tank from the top half of the flosser by rotating it anti-clockwise and filling in the required amount of water before properly closing it back again.

HOW TO USE:

- 1 To charge the device, check to ensure that the device and your hands are dry, and then connect the USB Plug to the 5V adapter.
- 2 After the device has ample amount of charge, you can disconnect it from the plug socket.
- 3 Fill the water tank with enough water for one use before selecting your desired mode by pressing the Mode button. You can also switch the modes while using the device.

The Product has three modes.

- a Cleaning: Deep clean and quick cleaning
- b Pulse: Gently massage the gums to promote blood circulation and improve gum health. It is suitable for sensitive teeth.
- c Surfing: Surge Point Spray. This mode is meant for massaging gums one by one, point-to-point, for a precise and an efficient rinse. For the first time users, using a dental appliance is suggested

Long Press the mode key for each gear to adjust to the pressure (30-120psi).

DEVICE CLEANING:

- 1 Open the Water tank cover and pour any remaining water out.
- 2 Press the button which flosses the teeth, this will drain the accumulated water.
- 3 Shut down the device before closing the water tank cover.
- 4 Dry the tooth cleaner with a dry cloth.

CARE INSTRUCTIONS:

- 1 Please use clean water of neutral detergent to clean the product. Using corrosive or abrasive detergent will damage the product.
- 2 Cleaning the Body:
 - 3 Clean the body with a dry cloth
 - a Do not Immerse the body in water for cleaning
 - b Cleaning the water tank: It can be cleaned with a dry or a wet cloth.
 - 4 Cleaning the nozzle: It can be cleaned with a dry or wet cloth, or removed and washed, then dried with a cloth.

DO'S AND DON'TS:

- 1 Do not use this product for purposes other than oral cleaning
- 2 Do not press the nozzle on the teeth or gums to avoid harming the gums.
- 3 Do place this product in a dry place.
- 4 Do not charge this product in hot or humid places.
- 5 Do not use water above a temperature of 40 degrees in the water tank or for cleaning.
- 6 Do not fill the water tank with mouth wash. Only use clean water.
- 7 Do consult a doctor before using the device in case you have had an oral surgery or suffer from any oral disease.
- 8 Do not bend, pull or twist the nozzle hose.
- 9 If not used for a long time, please unplug the device to avoid any risk.
- 10 This product is NOT suitable for children or infants.

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from the date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log onto to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

1. This warranty is void if:

- a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.