




Oral irrigator

Instruction manual

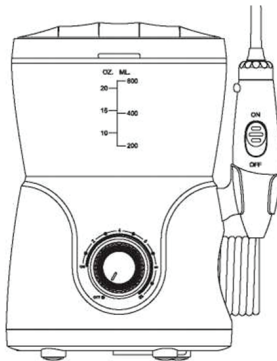

COMPLETE
ORAL CARE



600 ML RESERVOIR
CAPACITY



360 DEGREE
ROTATING TIPS



3 MINUTES SMART
TIMER

1 YEAR
WARRANTY

Thank you for choosing our oral irrigator!

You are on your way to better gum health and a healthier smile! With the combination of water pressure and pulsation, the Oral Irrigator removes harmful bacteria deep between teeth and below gumline where brush can't reach.


Traditional Clean

Traditional brushing can't easily reach hidden food debris and harmful bacteria.



Oral Irrigator Clean

Removes trapped food debris and harmful bacteria with high frequency pulse water floss. Flossing has limited access below the gumline

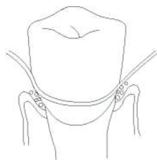

ACCESSORIES

Main Unit (Including Reservoir)	x 1
Instruction Manual	x 1

Tips

Ordinary High-Pressure Tip	x 3
Periodontal Tip	x 1
Orthodontic Tip	x 1
Dental Plaque Tip	x 1
Tongue Cleaner	X1
Hydraulic rotating toothbrush	X1

IMPORTANT SAFEGUARDS


When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

To reduce the risk of electrocution:

1. Always unplug this product immediately after using.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop in water or any other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately. To reduce the risk of burns, electrocution, fire, or injury to persons.
6. The product should never be left unattended when plugged in.
7. Close supervision is necessary when this product is used by, on, or near children.
8. Use this product only for its intended use as described in this manual. Do not use any other attachments.
9. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Please get in touch with our customer care.
10. Keep the cord away from heated surfaces.
11. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and dust.

12. Never use while sleeping or drowsy.
13. Never drop or insert any object into any opening or hose.
14. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
15. Connect this product to a properly grounded outlet only.
See Grounding Instructions.
16. Unplug this product before filling. Fill with water only. Do not overfill.



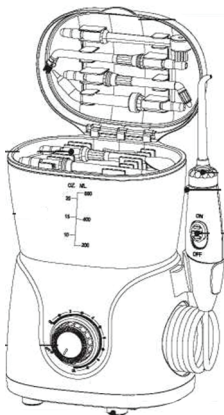



PRODUCT STRUCTURE



Ordinary High-Pressure Tip: 3 color ring. For hygienic purpose. Recommend using by 3 people respectively (family pack).


The other 4 special tips

- Periodontal Tip x 1
- Orthodontic Tip x 1
- Dental Plaque Tip x 1
- Tongue Cleaner x 1
- Hydraulic rotating toothbrush x 1

					
Ordinary High-Pressure Tip	Periodontal Tip	Orthodontic Tip	Dental Plaque Tip	Tongue Cleaner	Hydraulic rotating toothbrush


Instructions

Plug into socket. If the socket is controlled by a wall switch, make sure the switch is turned on.

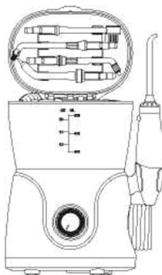


Inserting and removing tips

Insert the tip into the center of the knob at top of the dental jet handle. The colored ring will flush with the end of the knob if the tip is correctly locked into place. To remove the tip from the handle, press the jet tip eject button and pull the tip from the handle.

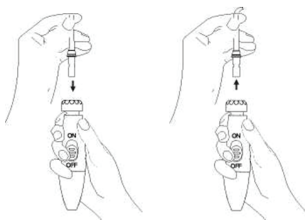


Adjusting The Pressure Settings

When you hear "click" sound, the switch has been opened. Then turn the knob clockwise to adjust the pressure from low to high. The machine will stop automatically after 3 minutes

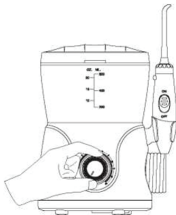

Using the Oral Irrigator

Lean low over the sink and place the jet tip in your mouth. Aim the tip towards your teeth. With your free hand, turn the on/off button on. Set pressure at desired setting. Direct the stream by turning the knob at the top of the handle.

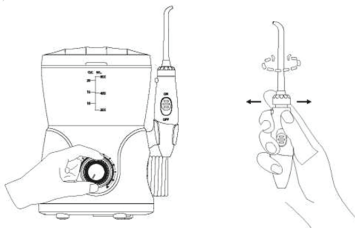


Recommended technique

Direct the jet stream at 90-degree angle to your gumline. Slightly close lips to avoid splashing but allow water to flow freely from mouth into the sink. For best results, start in the molar area (back teeth) working towards the front teeth. Glide jet tip along the gumline and pause briefly between teeth. Continue the process until all areas around and between the teeth have been cleaned.



Pause Control

You can temporarily stop the flow at any time by sliding the button downwards on the handle.

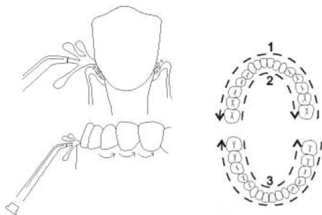

When finished

Turn the unit off. Dispose any liquids left in the reservoir. Unplug the unit.



INTRODUCTION TO THE JET TIPS

Periodontal Tip

The periodontal tip can direct water or liquid oral medicines in specific locations beneath the gum lining.



Ordinary High-Pressure Tip


Flosses deep locations within the oral cavity with a 360 degree range, supplemental^ covering places which are hard to clean, food scraps and harmful bacteria located deep within dental cracks or beneath the gum.

Tongue Cleaner Tip

Place the tongue scraper on the tongue before switching on the Oral Irrigator. Lightly place the tongue scraper on the centre-rear part of the tongue, lightly pulling it forwards while keeping it in contact with the surface of the tongue.

Dental Plaque Tip

Dental plaque is a group of bacteria that stick to the surface of the teeth yet has not mineralized. The jet tip is very effective at preventing dental plaque. The special-use dental plaque tip is for cleaning areas which is difficult to reach in braces and other dental devices.

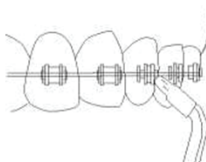

Orthodontic Tip

The special Orthodontic Tip is uniquely designed to simultaneously brush and rinse hard-to-reach areas around braces and other dental work.

To use the Orthodontic Tip, lean low over the sink and use the brush bristles to clean around dental work. Water will flow through the brush to rinse away any debris that is removed.

Hydraulic rotating toothbrush

The brush head is designed according to the dental scaling tools, which is different from the traditional toothbrush and just covers the whole tooth, thoroughly cleaning the surface of each tooth, the plaque hidden inside the tooth and the gingival crevice, and will not hurt And gums.



Cleaning

After using any special solution in the reservoir, rinse the unit properly.

Partially fill the reservoir with warm tap water, aim the jet tip into your sink, turn on the unit and let it run until the reservoir is empty. This will help keep the tip from clogging and let your unit to last longer.

Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser.

Before cleaning, please unplug from the electrical outlet.

Before exposing the dental jet to freezing temperatures, remove the reservoir and run the unit until completely empty.

Service Maintenance

1. After using our oral irrigator, be sure to unplug the power cord, and keep the tank upside down on a flat surface. When the device is completely dry, put it into a cool, dry place. It should be placed in the box or a similar storage when not in use for a long time.
2. If the machine malfunctions, please contact our customer care
3. Do not disassemble or repair the product on your own.
4. Contact our customer care if you have any problems with the product.

COMMON PROBLEMS AND SOLUTIONS

Problem	Cause	Solution
Jet Tips and Tongue Cleaner Leakage Between the top and handle knob	Tip not fully engaged	Remove tip and reinstall
Periodontal Tip Tear in soft tip	Tip is worn	Replace with new tip
Oral Irrigator Base hose to the handle leaks	Damage to the hose	Can be replaced by the seller for a self-service repair kit
Reservoir leaks	Missing water tank valve or water tank valve is upside down	Replace Reservoir valve (contact the seller); Turn the valve right side up
Inadequate pressure	Reservoir not positioned properly	Be sure outlet is working, test with another appliance (check electrical reset button and or wall switch)
Unit without water flow out	Motor not running	Fill out reservoir

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log onto to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised service.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customer@lifelongindia.com
www.lifelongindiaonline.com