



Rechargeable Water Flosser

Instruction



230 ML
WATER TANK



360 DEGREE
ROTATION



3 MINUTES AUTO
SHUT OFF

Safety Instruction

⚠ Dangerous

Prohibition: Do not block the head during use.
Prohibition: Do not use more than 10 minutes at a time.

1. Please read the manual carefully before use. The improper operation can cause injury or damage.
2. Do not disassemble this product. The product can get damaged or cause an electric shock.
3. Water poured into the tank should not be more than 40 degree.
4. If the product falls into the water when charging, please unplug the Power plug first.
5. If the power cord or plug is damaged, please replace the power adapter then use.
6. Do not put the product into the fire

⚠ Warning

1. Do not use the product for other purposes other than oral cleaning.
2. Do not use the nozzle pressure on the gums or teeth, so as not to damage the gums.
3. Do not under the product into the water to avoid the product failure.
4. Do not use this product in high temperature and humid environment
5. Do not continue to use the power cord if the power cord is damaged.


⚠ Note

1. Do not use this product without installing the nozzle.
2. Some mouthwash may have a certain corrosive. Do not put into the tank.
3. If a child wants to use this product, ensure that there is some supervision or guidance on the use.
4. Please keep children and babies out of reach, do not use as toys.
5. Patients with periodontal disease or have got an oral surgery in the last 2 months or are suffering from dental problems. Please consult a doctor before using


Advantage

- If used for more than 6seconds in same speed, the product locks the mode and remembers the mode when you restart it again
- Device should be turned off when charging.

Accessories

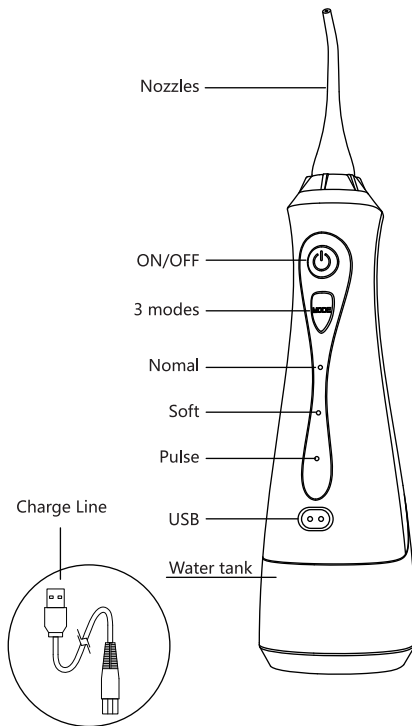


Structure of the device

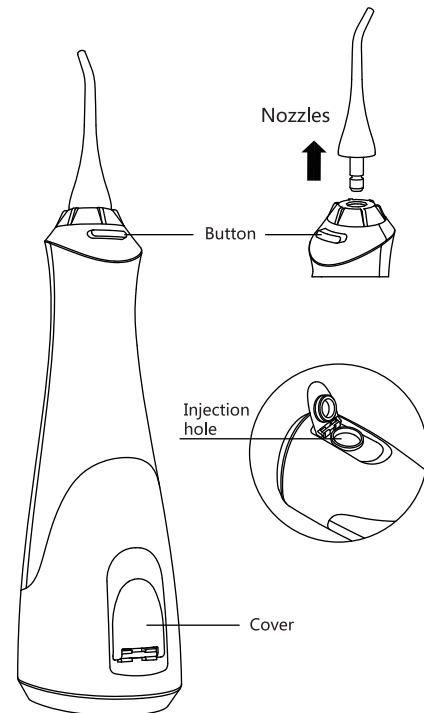

Use steps


1. Add Water

Open injection cover and set the product horizontally. Put under the faucet directly or use a glass to pour water in the tank.

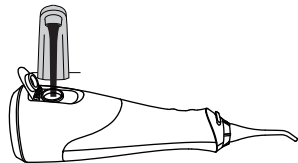

2. Clean Teeth

When cleaning the teeth, extend the nozzle into the mouth, aim at the teeth or gums, and open the mouth slightly so that the water flows smoothly.

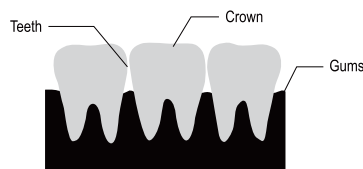

Daily Maintenance

- If you need to clean this product, please use dry cloth or wet cloth to scrub the body and nozzle.
- Do not immerse the body in the water.
- Do not use more than 40 degrees of hot water.
- Do not bend, stretch and twist the nozzle hose.
- It is recommended that nozzles be replaced every 6 months. Please dispose the water in the tank and place it in a dry and shady place.

FAQ

Q: Is it normal to use products that cause gums to bleed?

A: Gum bleeding is caused by inflammation caused by bacterial propagation between the teeth. Cleaning teeth each morning and night can help most users to slow down inflammation and bleeding.

If bleeding continues, please consult your dentist

Q: How to judge the effect of the tooth punch?

A: The role of the punch is to remove dental plaque between the teeth, helping to improve the health of the gums. Dental plaque is not visible to the naked eye, so if the following symptoms appear, that means that the product has been effective: overtime, the gingival bleeding situation is decreasing, the gingival tissue appears pink instead of red, the gums between the teeth appear to be sharp, rather than an inflamed arc-shaped. Usually after 24 weeks of use of the flosser, it will significantly improve the health of your gums. You can consult your dentist about your gum improvement.

Q: Can the teeth be cleared of all the things between the teeth?

A: As part of the daily oral care program, the effect of the punch is to clear the plaque between the teeth, which should remove most of the food **residue**. However, it is not guaranteed to be used to remove large, and sticky food residues.

Product specification

Net weight: 295g

Charge time: 2-4h

Power: 5W

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Tank capacity: 230ml

Input: 5V, 1A

Battery capacity: 1800mAh

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log onto to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :
Lifelong Online Retail Private Limited
For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com