




1 YEAR
WARRANTY

Smart Body Fat Scale

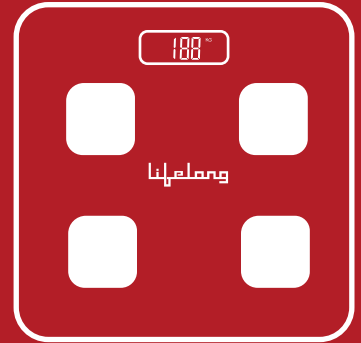

Manufactured for & Marketed by :
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com



150 kg



Health
Check



Multi User
App



Auto
On/Off



Battery Level
Indicator

INTRODUCTION:

Welcome to the Lifelong family.

Dear Customer,

Congratulations on purchasing Lifelong Smart Body Fat Scale and a warm welcome to the Lifelong family.

At Lifelong, we take pride in developing quality products for your home, health and happiness. Through constant innovation in our products, we have been able to develop a wide spectrum of products across Home, Kitchen, Grooming and Wellness categories. We strive to enrich the lives of the modern consumer every day. Keeping that in mind, we have offered all our products at a fair and honest price.

We hope that the Smart Body Fat Scale serves you and your family's needs without any hassles.

Thanks again for becoming a part of the Lifelong family. Enjoy using the product!

Warm Regards,
Team Lifelong

INSTALLATION:

- Insert two AAA Alkaline batteries into the battery compartment.
- Place the body scale on a firm, level and dry surface.
- Avoid places where potential water splashes may occur.
- It is advised to keep the scale away from daylight for better visibility of the display.




INSTALLING THE APP:

1. Please scan the QR code below to download the "Lifelong Smart" Home App.




2. If you don't have an app account, you can choose to register using an email address.
3. Go to the registration page, then click on "Get Verification Code" and enter it.
4. Set your password and click "OK" to register.

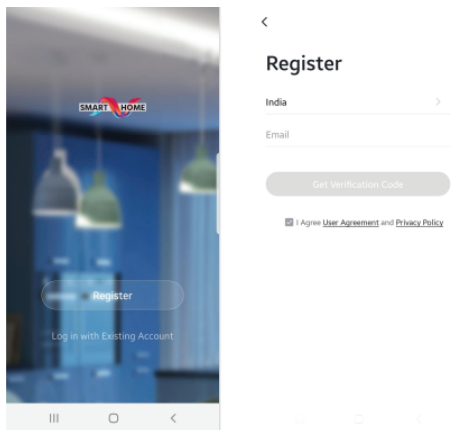


ADDING THE DEVICE:

1. Make sure your Bluetooth and Wifi are enabled on your mobile phone.
2. In the Lifelong Smart App home page, tap the icon "+" on the top right corner or click the "Add Device" button.

3. Select "Auto Scan" from the top carousel.
4. Press your foot on the body scale to activate it.
5. The body scale will be detected automatically in a few seconds and added.

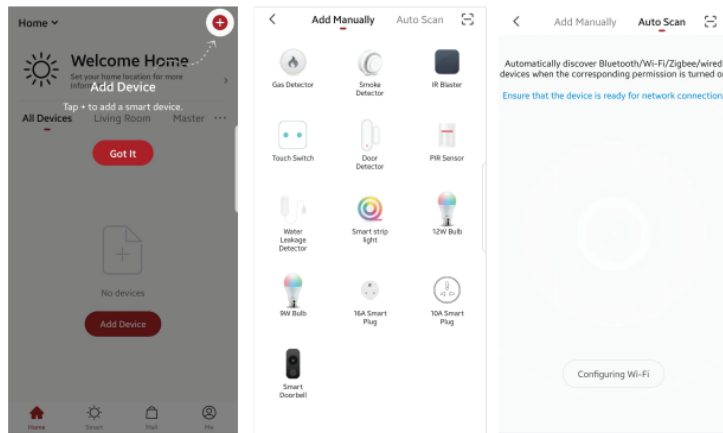

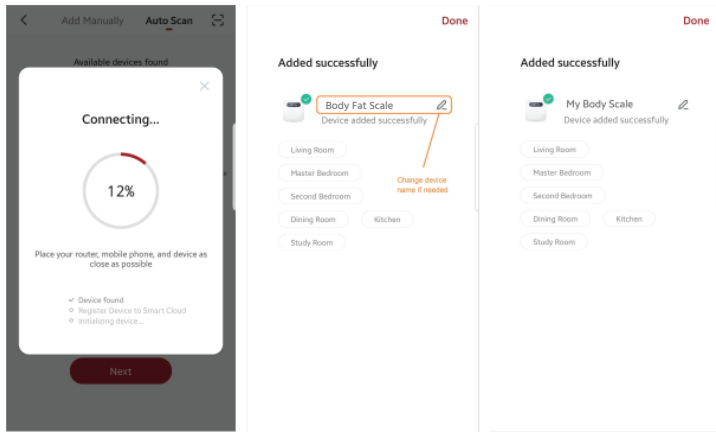

OPERATION:

1. You may now rename the device as per your choice.
2. You will then be prompted to select the unit of weight measurement (lb or kg).
3. In the subsequent steps, you are asked to enter your gender, date of birth and height along with a username.
4. To add more users, this step can be repeated from the Settings> Members> Add member.

5. On returning to the main screen of the body scale, you can now select the User and start measurements.

6. You can also set your target weight and target fat % on the main screen.

7. Once the measurement is completed, the data is transferred from the body scale to the mobile app. Please remember to keep Bluetooth enabled on your mobile during measurements.

8. The measured weight and body parameters can be viewed by scrolling down the main screen.

9. You can also view your measurements Trend and History from the bottom carousel.

10. History and Trend is stored for each user separately.

11. The battery level can be monitored from the app settings.


Note:

A. It is important to select the User before the start of the measurement. If not, the data of the other person will be recorded under the historic data of the user selected in the app. It is not possible to delete individual historic data entries.


B. Body weight can be measured even if the body scale is not connected to the mobile app.

But later when the body scale is connected to the app, all the data stored during the disconnected time will show up as “Unallocated Data” in the app. This data has to be deleted manually before starting a new measurement.


CHANGING/UPDATING THE SETTINGS:

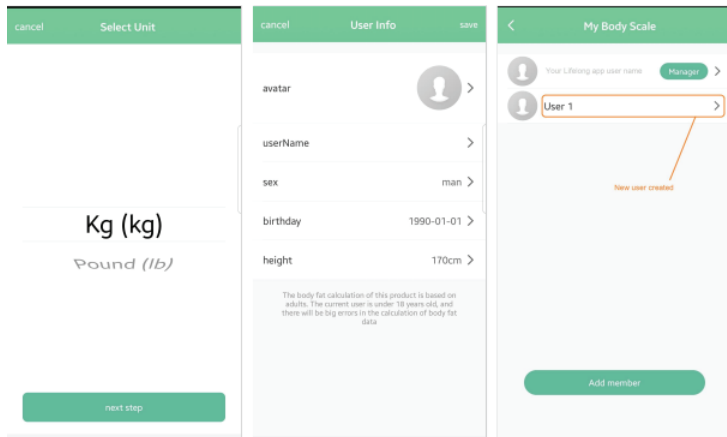

THE DETAILS ARE DISPLAYED ON THE APP ONCE MEASUREMENT IS DONE:

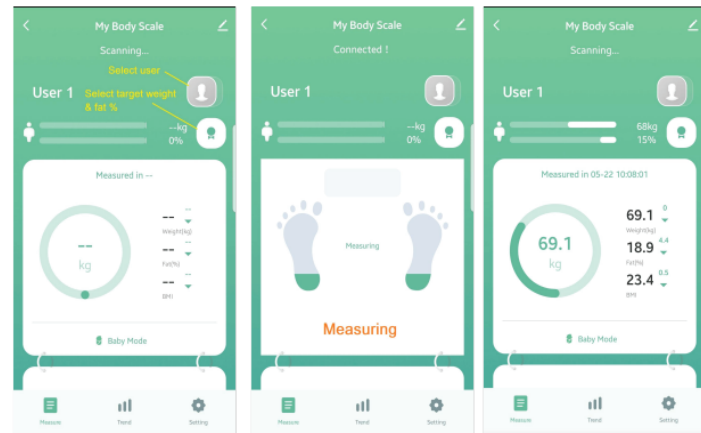


CHECK TREND AND HISTORY FOR USERS:

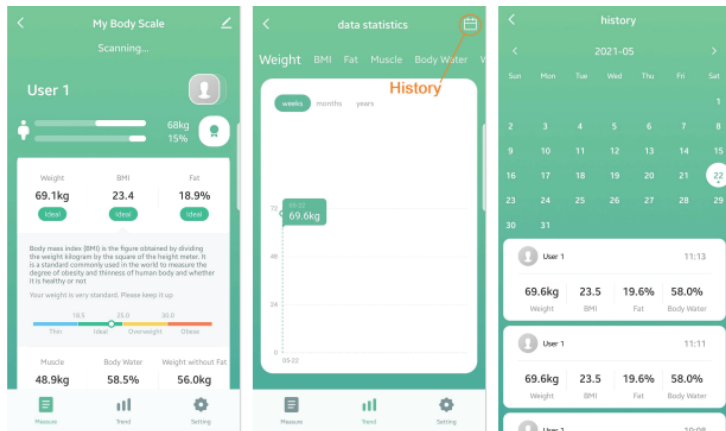

UNIT:

If users want to change units, please click <Set up>--<Unit switching> to select weight


unit :kg/lb/st

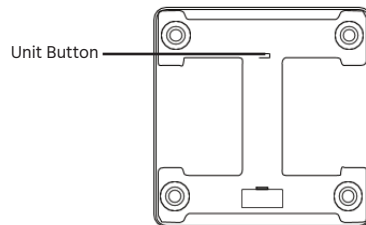


ERROR ICON OF THE SCALE:

Overload Or Measure In Error Operations:

the LED will indicate "Err", when scale is over max capacity. Please remove the weight to avoid any damage.

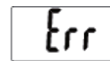

Low Battery : The LED will indicate "Lo", Please open the battery cover and replace the new battery with correct way.

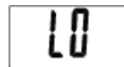

Faulty Measurement: Measurement fail indicates“Err1”, for below 2 conditions:

1. Body fat percentage is less 5% or over 50%, 2. Unsuccessful test.

SPECIFICATIONS:

- Power Source: 3V DC (Battery operated)
- Battery: 2 x AAA Alkaline Batteries
- Maximum weight capacity: 150kg
- Display: Liquid Crystal Display
- Backlight: Blue
- Working humidity: <85%
- Operating temperature: -10°C to +45°C
- Product dimensions: 280x280x24mm
- Package dimensions: 300x320x40mm
- Net weight -1400g
- Gross weight - 1700g

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.