

وسملانا

## DETACHABLE BODY MASSAGER LLM54



Powerful precision  
tapping massager



Removable and long ring  
handle for some hard-  
to-reach parts



2 removable massage  
head attachments

Dear Customer

Congratulations on your Lifelong Purchase. Our name stands for high-quality and thoroughly tested products. Please read these instructions carefully before use.

## HOW TO USE

The Lifelong Detachable Body Massager can give you an effective full body massage in the comfort of your own home. The sturdy design and detachable feature ensure you can reach hard to reach areas of the body as well with ease.

Step 1: Ensure that the device is switched off before you plug in.

Step 2: Choose between 2 different massage head options – Multiple point and rounded depending on your personal preference

Step 3: Detach the massager with the two buttons on each side of the massager if you want to use it on easily accessible areas or use the full version for areas such as back etc.

Step 4: Put your hand under the strap for a firm grip

Step 5: Switch on the device by sliding the switch to ON position

Step 6: The intensity and the speed of the massager can be adjusted with the button and the slider given on the upper left and right side of the device. You will also see LED changing intensity from high to low as you adjust the device according to your required comfort.

Step 7: Guide the device to the area of the body to be massaged. You can stand, sit or lie down and use the massager

## FUNCTION

1. Split structure design, light and handy, great for a massage anytime, anywhere.
2. The main part of the massager can be used alone, relying on the locking strap fixed in the palm to make the massage simpler.
3. Massager's handle is convertible. When attached one can reach hard to reach areas like back etc. and when detached it becomes much smaller for general use.
4. Continuously variable speed regulation can effectively stimulate the nerves and muscles for a refreshing massage.
5. Point massage helps relieve muscle pain.

## SPECIFICATIONS

Input voltage: AC220-240V

Frequency: 50-60Hz

Rated Power: 25W

## Consumer Warranty Card

Dear Customer,  
Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 6 months

Please note:

Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log onto to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.


Step 6: Use on a particular area only for 15 minutes and then switch to another. Massager use is recommended to be upto 30 minutes to avoid overstimulation of nerves.

Step 7: Unplug the device and let it cool down after use. Keep the device dry.

## FEATURES

### MASSAGE HEAD

Two different shapes of the massage cover design can help concentrate on specific areas.

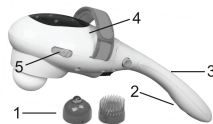

### THREE COMBINATIONS

Three different combinations mode, help switch as per requirement of the user. The locking strap design helps the product stay in place while massaging.



### MASSAGE STRENGTH ADJUSTMENT

Choose between the different massage strengths from soft to harder.

### ONE KEY SWITCH DESIGN

The switch has three gears, the first gear to start the vibration and infrared, and the third gear to turn off the internal power supply.

## Terms and Conditions

### 1. This warranty is void if:

- a. The completely filled warranty card is not presented at the time of servicing the product.
- b. The product is not operated according to the instructions given under the user manual.
- c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
- d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
- e. Product label specifying the model number, serial number and production code has been removed and altered.
- f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
- g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
- h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.

### 2. Repairs or replacements will be carried out by authorized service provider.

### 3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.

### 4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply notwithstanding the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

Made in China

For queries and complaints please contact:

customercare@lifelongindia.com

www.lifelongindiaonline.com