


Lifelong

Lifelong
Rechargeable
Cushion
Massager



Infrared
Heat Therapy

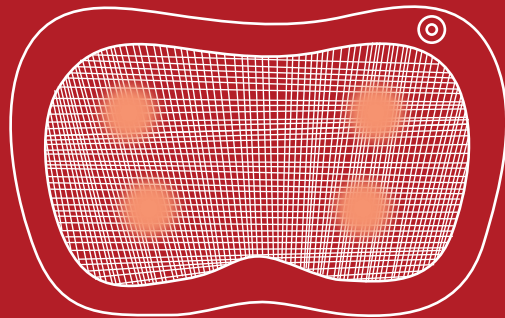

Full Body
Relaxation



Shiatsu Body
Massage



Rechargeable

1 YEAR
WARRANTY

INTRODUCTION:

Welcome to the Lifelong family,

Dear Customer,

Congratulations on purchasing Lifelong Rechargeable Cushion Massager and a warm welcome to the Lifelong family.

At Lifelong, we take pride in developing quality products for your home, health and happiness. Through constant innovation in our products, we have been able to develop a wide spectrum of products across Home, Kitchen, Grooming and Wellness categories. We strive to enrich the lives of the modern consumer every day. Keeping that in mind, we have offered all our products at a fair and honest price.

We hope that the Cushion Massager serves you and your family's needs without any hassles.

Thanks again for becoming a part of the Lifelong family. Enjoy using the product!

Warm Regards,

Team Lifelong

WARNING:

If the unit is not working or not functioning properly, you should immediately stop using it.

Do not use it for a long time. Not more than 30 minutes each time

When the power plug is ON, do not pull the cable with bare hands

Do not place it in high humidity or use in the bathroom.

Do not immerse in water or bring the power supply near water, so as to not cause a short circuit or electrical shock.

Children should use it under the care of adults.

Inappropriate use of this product

- During menstrual period and pregnancy women, patients with severe heart disease, Hemorrhagic disease, sepsis, hives and other patients.
- Pay attention to Children and mentally handicapped people with disabilities should use this under supervision.
- The machine with heating function, is not sensitive to heat
- Please don't put too much pressure or stand on the phone, to slow down or stop it, as it can lead to motor damage of the machine.
- Please do not remove the protective cloth.
- In case of damage or inflammation of the skin, do not use.
- Please drink 500 ml water or tea after massage, so as to promote the body metabolism.
- If it is found that the machine shell is broken or the cloth is torn, do not insert your fingers into rotation massage head bottom, in order to avoid the task.

Use environment

- Do not use the product in humid or dusty environment to avoid electric shock or machine failure.

- Do not place the product near the stove or in direct sunlight.

*Note: Do not use it during lightning strike.

The battery safety

- Do not connect the metal conductor with the battery poles or connect the battery ends to avoid short circuit of the battery and body injuries such as burns caused by battery overheating.

- Do not expose the battery to high temperature or to heating equipment, such as sunshine, heater, microwave, oven or water heater, etc. Overheating battery may cause explosion.

- Do not disassemble or modify the battery, insert a foreign body, or immerse the battery in water or other liquids, so as to avoid battery leakage, overheating, fire or explosion.

- If the battery leaks, please do not let the skin or eyes come into contact with the leaking liquid. If it comes into contact with the skin or eyes, please rinse with water immediately and ask the doctor for medical treatment.

- If the battery is discolored, deformed, abnormal heating and other abnormal phenomena during use, charging or storage, please stop using and replace the battery: Do not throw the battery into the fire, otherwise it will cause the battery to catch fire and explode.


- Please dispose of the battery according to local regulations, do not dispose of the battery as domestic garbage, improper disposal of the battery may lead to battery explosion

- Do not drop, squeeze or puncture the battery, so as to avoid the battery being subjected

to large external pressure, which may lead to short circuit and excessive heat inside the battery.


- when the product use time is significantly shorter than normal time, please replace the battery.

PRODUCT APPEARANCE STRUCTURE:

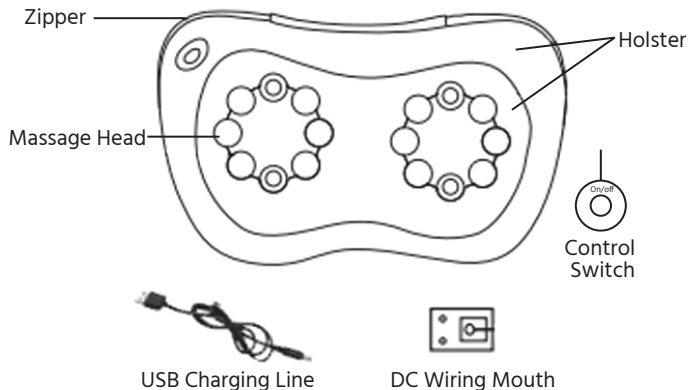

CHARGING:

1. Before first use, charge for about four hours.
2. In the process of charging, charging indicator lights will flash.
3. Once charging is complete, the charging light is normally on.
4. When using the product, if the battery is low, charging indicator lights will flash.

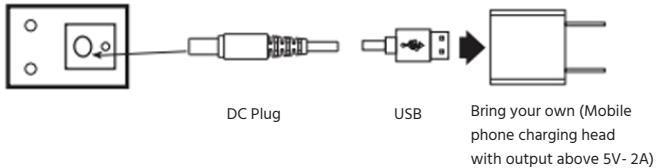

OPERATING INSTRUCTIONS:

Home use: connect the mains plug and use.

Car use: the power adapter connects to the vehicle power supply.

- Press Power button once on the top left corner to start the massage and heating function simultaneously
- Pressing the button again to start kneading massage function which will automatically change every minute to massage
- Press the button the third time, the product restarts the heating function again, like a new cycle.

- Press the button for 3 seconds to turn off the product the system returns to standby.

Note: The system features automatic timing 15 minutes, we suggest using the time off every half hour.

FEATURES:

This product uses a unique design:

1. Four heat kneading massage heads that massage the deeply fatigued muscle tissue, relieve pain, relax, allowing you to more fully and comfortably enjoy massage.
2. A controllable infrared heating function that can produce heat effect.
3. Automatic reversing feature that allows you to experience two different massage
4. Automatic 20 minutes time function.

CAR & HOME MASSAGE PILLOW:

Using high-tech computer digital technology, humanoid hand massage edition which combines advanced ergonomic theory, accelerate blood circulation, rapid elimination of fatigue, neck pain, stiffness, rapid lifting of job pressure and restores a happy mood.

CLEANING AND MAINTENANCE:

- Before cleaning, pull the plug from the socket.
- Please put the machine in dry area to prevent dust from damaging the product.
- Avoid contact with water or other liquid to the machine.

- Please do not use volatile oil, bleach or other corrosive solvents to wipe the machine.
- This machine is specially design and manufactured so that no maintenance is required, also need not add lubricating oil regularly.

TECHNICAL DETAILS:

Rated voltage	AC100-240V/DC12V
Rated frequency	50-60HZ
Rated power	25W
Rated time	15mins

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

- This warranty is void if:
 - The completely filled warranty card is not presented at the time of servicing the product.
 - The product is not operated according to the instructions given under the user manual.
 - Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - Product label specifying the model number, serial number and production code has been removed and altered.
 - Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
- Repairs or replacements will be carried out by authorized service provider.
- During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised service provider.
- All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
- Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
- Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
- However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not withstanding the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
- No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com