

Manufactured for & Marketed by:

Lifelong Online Retail Private Limited

For queries and complaints:

please contact: customercare@lifelongindia.com

www.lifelongindiaonline.com

Toll Free +91-9667499699



Electric Handheld Full Body Massager



USER MANUAL

INTRODUCTION:

Dear Customer, Congratulations!

We welcome you to the Lifelong family, as you welcome our Electric Handheld Full Body Massager to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Electric Handheld Full Body Massager.

We know we can't!

Looking forward to building our bond for life, cheers!

Thank you,

Team Lifelong


TABLE OF CONTENTS :

1. About the product
2. Product parts
3. Safety precautions
4. Using the product
5. Maintenance
6. Terms & Conditions
7. Warranty Card

ABOUT THE PRODUCT:

The Lifelong Electric Handheld Full Body Massager is a rejuvenating device that gives you great value for money. The massager kneads and relaxes your muscles, letting the days worries melt away.

PRODUCT PARTS:

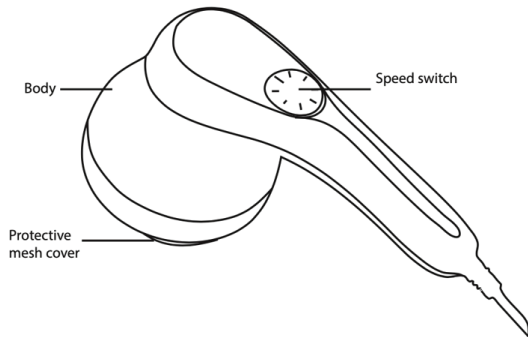

SAFETY INSTRUCTIONS:

- Use this appliance only for its intended use as described in this manual.
- Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Do not carry this appliance by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair.
- Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- Do not reach for the appliance that has fallen into water. Unplug immediately.
- Do not use while bathing or while showering.
- Do not place or store appliance where it can fall or be pulled into a tub or sink.
- Do not place in water or other liquids.
- If the supply cord is damaged, you must get it repaired through our customer care.
- This product should not be used by pregnant women, or people with skin disease, inflammation of the skin, heart disease and other disease.
- Children and disabled persons should not use the machine.
- Do not use for more than 15 minutes in the same position.
- It is suggested to not use the massager within an hour of food intake.
- Please make sure the machine is off before changing the massage head.

FEATURES:

- Heavy duty motor: Extremely powerful energy efficient copper wired motor encased by premium quality abs for superior performance
- Dual operation: Perfect choice as body massager & toning for weight loss. Relax your aching muscles while targeting those extra fat pockets to help reduce weight and improve stamina at the same time
- 4 massage heads: Targets different parts of body for customized massage
 - a) Scraping: scraping effect
 - b) wavy: acupuncture c) ball: use with massage oil d) micro fibre: dry skin removal from feet
- Additional mesh cover: The protective soft mesh fabric head is provided only to avoid hair entanglement during massage. Do not use with micro-fibre head as fabric is delicate & easily prone to damage
- Reliable features: Extra-long 1.6M cord for maximum portability; ergonomically designed and gripped handle for easy movement; all-in-one step less speed regulator knob to set speed between low to advanced.

HOW TO USE:

1. Put the plug into the socket.
2. Turn on the speed switch to adjust the different massage strengths. Rotate the speed control switch clockwise to gradually increase the massage strength. Rotate anticlockwise to decrease the massage strength

3. After using, turn the switch "OFF", and then pull out the plug.

MAINTENANCE:

1. After using, place the product in a clean place, with no direct sunlight.
2. Do not pull the cord or pull the product with cord.
3. Make sure that the product does not come in contact with water.
4. Please use a soft cloth to clean the machine.

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

Terms and Conditions:

1. This warranty is void if:

- a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
 3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
 4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.