

Manufactured Market by:
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com

Lifelong

Foot and Leg Massager



Kneading, rolling and
vibrating actions



Preset and
manual modes



Washable
foot cloth

**1 YEAR
WARRANTY**

INTRODUCTION:

Dear Customer, Congratulations!

We welcome you to the Lifelong family, as you welcome our Foot And Leg Massager to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Foot And Leg Massager.
We know we can't!

Looking forward to building our bond for life, cheers!

Thank you, Team Lifelong


TABLE OF CONTENTS :

1. About the product
2. Marking of product & Additional Features
3. How to use
4. Device Cleaning
5. Care instructions
6. Cord Safety
7. Safety for usage/ Dos and Don'ts
8. Technical data

ABOUT THE PRODUCT:

The Leg Massager is a calming device that gives you great value for money. The massager kneads and relaxes your muscles, letting the day's worries melt away. With more than 100 innovative products available, it's no wonder that millions of customers continue to choose Lifelong! It has been strongly trusted by loyal users for over five years. Made specially for you and your relaxation.

PRODUCT INFORMATION:

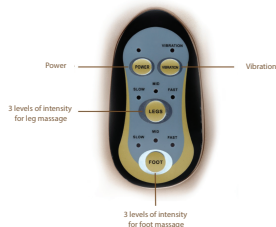

Lifelong Foot Massager is based on customizable massage with a 3-modes of Vibration and Intensity based reflexology massage therapy. Helps in relieving muscle tension and aids better blood circulation. 15-minute Auto Shut-off. The Lifelong Leg and Foot Massager is aptly designed to provide electronic simulation on the required pressure points at the sides of the sole and helps you relax, Maintains Blood Pressure, strengthens muscles, improve balance and flexibility.

HOW TO USE:

The Lifelong Foot and Calf Massager can give an effective foot massage to your feet, ankles and calves all at the same time.

STEP 1: Ensure that the device is switched off before plug in.

STEP 2: Put your feet and calves in the black cloth area.

STEP 3: Switch on the device by pressing the ON button until the light above it turns on.

STEP 4: You can change the settings by clicking on the buttons provided.

STEP 5: The device has 3 buttons to choose from, among Vibration, Foot and Legs.

STEP 6: Press the mode you prefer once, twice or thrice based on your preferred strength.

STEP 7: To stop the massage, press the mode key again.

DEVICE CLEANING:

To clean the unplugged device after it has cooled down, wipe with a slightly damp cloth, and ensure that water doesn't get inside the device.

The foot cloth should be removed from the zipper and washed separately.

- Never operate the device especially if there are children in the vicinity.
- Do not use with animals or pets.
- Do not cover the massager with the damp cloth when it is in an operating condition, it may result to a subsequent breakdown or fire.
- Do not sleep while using the massager.
- Never plug the device into the mains power supply if the cable is damaged.
- Never hold or carry the device by mains cable.

TECHNICAL DATA/SPECIFICATIONS:

Features	3 Customized Massage Options
No. of Motors	4
Size (cm)	54 x 53 x 49
Net Weight	10.5 kg
Gross Weight	13.5 kg
Voltage	220-240V AC
Frequency	50 Hz
Wattage	80W

Consumer Warranty Card

Dear Customer,
Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log onto to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

CARE INSTRUCTIONS:

1. Ensure that the device is switched off before you plug in.
2. Connect the device only to the type of voltage stated on the nameplate.
3. As with every electrical device, this massage device should also be used carefully and cautiously in order to avoid the risk of electric shock
4. Unplug the mains cable immediately after use.
5. Prevent any contact with water (except when cleaning with a slightly damp cloth).
6. If water gets inside the casing, unplug the device immediately and contact customer service.
7. The device and in particular the mains cable should be checked for any damage before each use.
8. If a fault should occur during use, switch off the device immediately and unplug it from the socket.
9. The device is intended only for the purpose described in these instructions for use. The manufacturer can accept no liability for injury arising from incorrect or careless use.

CORD SAFETY:

1. Check that the voltage corresponds with the specification indicated for the massager.
2. Do not disconnect the AC plug from the electrical outlet roughly.
3. Do not connect or disconnect the AC plug from the electrical outlet with wet hands. It may result in electrical shock or machine trouble.
4. When disconnecting the AC plug, pull it out by holding the plug not the cord.
5. Do not place the AC cord under the massager or any heavy object.
6. Do not wind the AC cord around the massager as it may damage the cord and result in fire or electrical shock.
7. Do not operate the massager with a damaged cord or tamper with the cord. Contact the service center for inspection.
8. Use only attachments recommended by us.
9. Do not use the massager if the electrical outlet is excessively loose.

DO'S AND DONT'S:

- Do not use the massager while you're in the washroom
- Do not use the massager on wet body, this may cause short circuit or electric shock
- Those suffering from cancer, osteoporosis, tumours, high fever or pregnancy must consult a doctor before use.
- Do not use the massager if it is in a bad condition

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not withstanding the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.