

# lifelong

## AIRBAG FOOT MASSAGER

User Manual



Built for Relaxation. Built for India.

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

For queries and complaints:  
please contact: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)

## INTRODUCTION:

Dear Customer,

Congratulations! We welcome you to the Lifelong family, as you welcome our Lifelong Airbag Foot Massager to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Lifelong Airbag Foot Massager. We know we can't!

Looking forward to building our bond for life, cheers!

Thank you,  
Team Lifelong



**LLM432**

**TABLE OF CONTENTS:**


- 1. About the Product
- 2. Product Marking
- 3. Description on Keys
- 4. Features
- 5. Instructions
- 6. Safety Precautions
- 7. Troubleshooting
- 8. Specifications

**ABOUT THE PRODUCT:**


The Lifelong Foot Massager is ergonomically designed for maximum comfort. It offers a deep-kneading Shiatsu massage, heated, squeeze massage. It can relieve chronic foot discomforts from plantar fasciitis, and help with neuropathic foot numbness.

**PRODUCT MARKING:**

**Product Features**

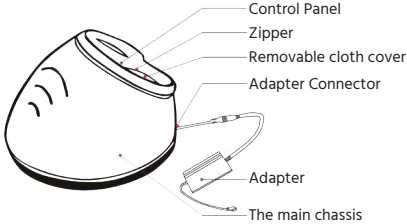

**Control Panel**

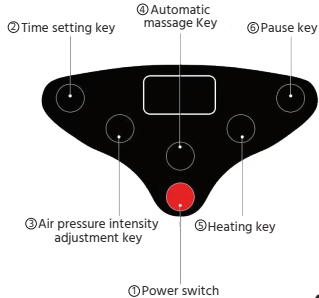

## DESCRIPTION OF KEYS:



POWER

### 1. Power switch

Turns the power on or off.

Default massage run time is 15 minutes, over 15 minutes will automatically turn off the power.



TIME

### 2 Time setting key

Turn on the power, you can set the massage time, press the time by 5 minutes, the longest massage time can be set to 30 minutes, the shortest is 5 minutes.



STRENGTH

### 3. Air pressure intensity adjustment key

Power on default is mid-range (L2).

Total weak (L1) - Mid-range (L2) - Strong (L3) Three-speed adjustable. The intensity is displayed in the display window (L1, L2, L3).



MODE

### MODE

4 6 massage programs can be adjusted

Built-in 6 automatic programs, respectively (P1, P2, P3, P4, P5, P6).

**P1:** For fitness massage function, (soles of the foot roller massage and air pressure massage at the same time work).

**P2:** For the health massage function; soles of the feet roller massage and air pressure massage at the same time work).

**P3:** For leisure massage function (foot roller massage and air pressure massage alternately work).

**P4:** For the foot roller independent massage (pinch feet).

**P5:** For the pressure massage (cover their feet).

**P6:** For the warm foot function (for heating the soles of the feet).

**Mode is displayed in the display window (P1, P2, P3, P4, P5, P6).**



HEATING

### 5. Heating construction

Turn on or off the heating function of the soles of the feet, the display window shows (01) that heating is turned on, the display shows (00) meaning heating is off.








PAUSE KEY

### 6. Pause key

Pause or turn on the running massage program.

## FEATURES:

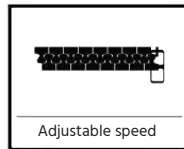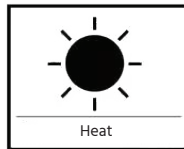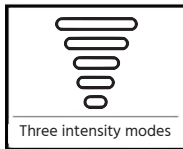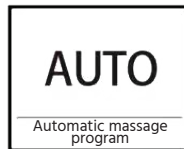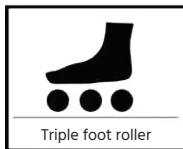

## INSTRUCTIONS:

1. Connect the product and the adapter, and then connect the power plug to the outlet.
2. Sitting in the chair for their height, feet into the foot machine inside.
3. Select the desired massage program.
4. Unplug the power cord after use.



## SAFETY PRECAUTIONS:

Read the safety information carefully before use and follow the instructions.

### NOTE

1. Recommended use time is 15 minutes, repeated use may lead to product overload.
2. Do not massage the wound or the injured body parts.
3. Users with sensitive skin should pay attention to not compress parts that may cause burns.
4. After each use, please turn off the power switch in a timely manner and remove the power plug.
5. Do not use the product near electric heaters or other heat sources.
6. Do not use excessive pressure on the product, as it may cause failure.
7. Keep the product away from humidity or high temperature.
8. Do not touch the product with wet hands or connect the power supply with wet hands.
9. Never pour water over the product, so as not to cause failure in case of dropping water, please immediately turn off the power switch and unplug the power plug from the outlet.
10. Do not dismantle the product.
11. Do not stand on the product.
12. Prevent children or users who may not be able to take care of themselves without authorization to operate the product.  
(Please ensure children and elderly people are accompanied by adults to use of products)
13. Please follow the doctor's advice if you are undergoing treatment or belong to the following cases:

- Cancer patient
- In vivo with metal devices, or with artificial limbs.
- Massage parts of wounded people.
- Chronic disease or patients with serious illness.
- Women having Menstrual period and pregnant women.
- A person with a heart or a mental disorder.
- A person who is being treated by a doctor.
- People with bone hyperplasia.
- A person with vertebral disease or bone fragility.

## TROUBLESHOOTING:

Fault type	Cause of issue	Maintenance methods
Unable to start the product	The power supply or adapter not connected	Replace the power supply or adapter
	Poor contract in the power plug	Re-insert the power plug, and then confirm
	Massage for more than 15 minutes	Turn the power switch back on
Suddenly stopped in use	Massage time is too long	Turn off the power, over 1 hour and then use
Use a faint sound	The motor carries the sound generally	Belong to the universal sound, non-failure
	Use of time is too long	

## SPECIFICATIONS:

Product type	Electric massage	Power	40W
Model name	LLM432	Suggested time	15 min
Voltage	AC 100-240V~,50/60Hz, 1.5 A	Weight	4.1 kg

## Terms and Conditions:

1. This warranty is void if:
  - a. The completely filled warranty card is not presented at the time of servicing the product.
  - b. The product is not operated according to the instructions given under the user manual.
  - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
  - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
  - e. Product label specifying the model number, serial number and production code has been removed and altered.
  - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
  - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
  - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

## Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

## Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)

## Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

## Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log on to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.