

Manufactured for & Marketed by :
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongindia.com
www.lifelongindiaonline.com

lifelong

LEG MASSAGER

User Manual



Built for Relaxation. Built for India.

INTRODUCTION:

Dear Customer,

Congratulations! We welcome you to the Lifelong family, as you welcome our Lifelong Leg Massager to be a part of your home.


At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Lifelong Leg Massager.

We know we can't!

Looking forward to building our bond for life, cheers!

Thank you, Team Lifelong

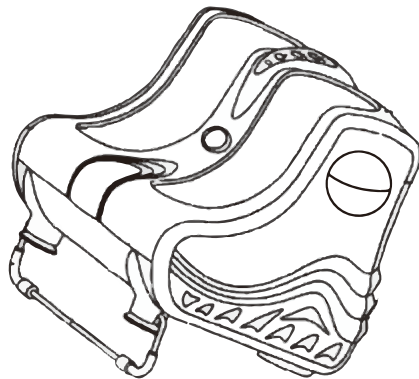

LLM747

TABLE OF CONTENTS:

1. About the Product
2. Marking of Product
3. Product Features
4. How to Use
5. Dos And Don'ts
6. Device Cleaning
7. Product Specifications
8. Troubleshoot Chart
9. Terms & Conditions
10. Warranty Card


ABOUT THE PRODUCT:

The Lifelong Leg Massager is a calming device that gives you great value for money.

The massager kneads and relaxes your muscles, letting the day's worries melt away.

Made specially for you and your relaxation.

MARKING OF PRODUCT:

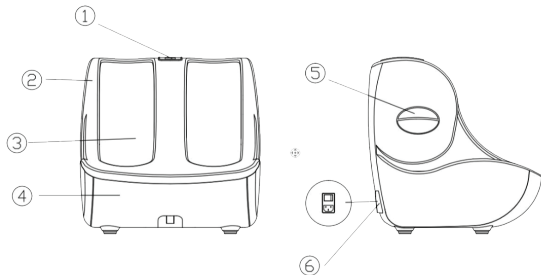

(1) Central surface: computer panel - control the massage function of the product

(2) Front upper and inner lining: kneading roller + massage plate - to achieve acupoint massage on the back of the leg and on both sides (3) front outer surface: cloth cover (with zipper) - detachable for cleaning

(4) Inside and outside the front bottom: kneading roller + massage plate + vibration - to achieve multi-functional acupuncture points in the foot and ankle joint area massage


(5) Both sides of the surface: handles - easy to move products

(6) Right side of the back: power switch - the power cord can be disassembled independently

PRODUCT FEATURES:

1. It can relax the muscles of the legs and feet.
2. It can stimulate acupuncture points and eliminate fatigue.
3. Can promote blood circulation.
4. It can relax tense nerves and reduce stress.
5. Can eliminate leg pain.
6. It can reshape the perfect curve of the legs.
7. Removable and washable cloth cover to keep you clean all the time.

HOW TO USE:

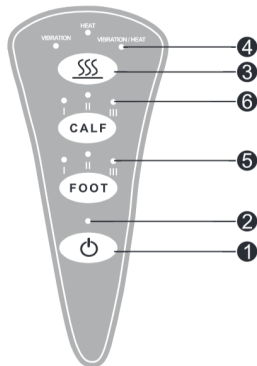

- 1.power button
- 2.indication light
- 3.vibration or heat button
- 4.indication light of vibration and heat
- 5.indication light of foot massager speed
- 6.indication light of leg massager speed

TO ACTIVATE THE KNEADING DISCS ONLY

1. Place the lower section of your legs into the fabric-covered slots. See Picture.
2. Press the "CALF" button and the default mode will be "I", press it again or twice, and mode will switch to "II" or "III". Different kneading modes for each level.
3. Press the "FOOT" button and the default mode will be "I", press it again or twice, and mode will switch to "II" or "III". Different kneading modes for each level.
4. Press the "Vibration/Heat" button to start function of vibration, press it again, the model will switch to heat, or press it twice, the model will switch to heat and vibration.
5. Press the "POWER" button to stop the kneading massage.



Open or close bracket according to the height of fee.



insert the adapter to the massager's back socket by correct direction, after plug connect power, the indication light of adapter will on.



stop use: please press the off button, then cut the power remove the adapter.



Choose chair be equal to stature, put feet into massage slot.



If the electricity was cut off in the course of using ,step one: drive up your knee.



If the electricity was cut off in the course of using ,step one: turn your feet side and move easily.

DOS AND DON'TS:

- 1. Ensure that the device is switched off before you plug in.
 - 2. Connect the device only to the type of voltage stated on the nameplate.
 - 3. As with every electrical device, this massage device should be used carefully and cautiously in order to avoid the risk of electric shock
 - 4. Unplug the mains cable immediately after use.
 - 5. Prevent any contact with water (except when cleaning with a slightly damp cloth).
 - 6. If water gets inside the casing, unplug the device immediately and contact customer service.
 - 7. The device and in particular the mains cable should be checked for any damage before each use.
 - 8. If a fault should occur during use, switch off the device immediately and unplug it from the socket.
 - 9. The device is intended only for the purpose described in these instructions for use.
- The manufacturer can accept no liability for injury arising from incorrect or careless use.

DEVICE CLEANING:

- 1.To clean the unplugged device after it has cooled down, wipe with a slightly damp cloth, and ensure that water doesn't get inside the device.
- 2.The foot cloth should be removed from the zipper and washed separately.

PRODUCT SPECIFICATIONS:

- VOLTAGE: AC220-240V~50Hz 1.6A
- POWER: 80W
- RATED WORKING TIME: 15 MINUTES
-

A/C ADAPTOR

- INPUT : AC220-240V~50Hz 1.6A
-

TROUBLESHOOT CHART:

If you encounter problems when using this product, please refer to the following tips and possible solutions. If the product is still in a faulty condition, please take the product to an authorized technician for inspection and repair.

Abnormal situation	Possible reason	Solution
Not working	1. Not powered 2. The function keys are not turned on	1. Plug in the power supply or check if the power supply is well connected 2. Turn on the on/off button
Sudden stop of operation	1. Power off 2. Massage for more than 15 minutes 3. If the continuous use time is too long, the product will be protected from overheating. 4. If the force is too large, the massage ball will be blocked.	1 Check whether the power interface is connected well 2. Turn the on/off button back on 3. Turn the on/off button back on after the product cools down 4. Reboot on/off button

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised service.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customer care@lifelongonline.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.