ومملعهنا

WIRED GUN MASSAGER

User Manual



Built for Relaxation. Built for India.

Manufactured for & Marketed by : Lifelong Online Retail Private Limited

For queries and complaints: please contact: customercare@lifelongindia.com www.lifelongindiaonline.com

INTRODUCTION:


Dear Customer,

We know we can't!

Congratulations! We welcome you to the Lifelong family, as you welcome our Lifelong Wired Massager Gun to be a part of your home.

At Lifelong, our products are carefully built for your needs, which is why we have crafted a diverse range of innovative products, across a spectrum of categories -Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living. We are sure you can't wait to start using your new Lifelong Wired Massager Gun.

Looking forward to building our bond for life, cheers! Thank you, Team Lifelong

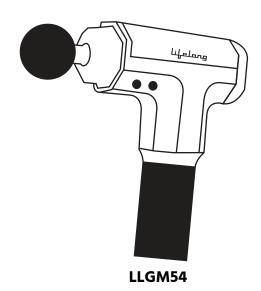

2

TABLE OF CONTENTS:

- 1. About the product
- 2. Safety Precautions
- 3. Operation Instructions
- 4. Changing the Massage Heads
- 5. Maintenance & Care
- 6. Dos and Don'ts
- 7. Terms & Conditions
- 8. Warranty Card

ABOUT THE PRODUCT:

The Lifelong Wired Massager Gun is a portable vibration device. It can effectively reduce the pain and relieve the discomfort in your muscles.

SAFETY PRECAUTIONS:

- · For adult use only.
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the Hyper volt only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Use only the applicators that provide the best desired effect.
- Bruising can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair or any other body part away from the shaft and backside of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the Massage Gun.
- Do not immerse in water or allow water to enter the ventilation holes of the Massage Gun.
- Do not drop or misuse the Massage Gun.
- Carefully examine device and battery before each use.
- Do not tamper or alter Massage Gun in any way.
- Never leave the Massage Gun operating or charging unattended.

OPERATING INSTRUCTIONS:

Operating:

- 1. Power device on by moving the power switch on the bottom of the battery to the ON position. Power device off by moving the power switch to the OFF position.
- 2. With power switch in ON position, press one time to start the device, press the 2nd time to change gear and press the 3rd time to switch it off.
- 3. With unit off, fully insert applicator by pressing in. Remove by pulling it out.
- 4. Massage desired body part at the desired speed applying pressure as needed.

Massage Heads:

The massager comes with 4 detachable attachments so that users can target all parts of their body.

- ROUNDHEAD: Suitable for massaging large muscle groups such as abdomen, glutes and legs. Great for muscle relaxation before and after exercise.
- FLATHEAD: Suitable for large muscle groups such as inside and outside muscles of the abdomen, chest and legs
- U-SHAPE HEAD: Suitable for shoulders, neck and bilateral muscles do not use directly on the cervical spine or spine.

 CONICAL HEAD: Suitable for applying direct pressure on joints, deep tissue, trigger points and release muscles knots.

Choose the speed and massage heads that best suit your needs to more precisely target different body parts and muscle groups.

CHANGING THE MASSAGE HEADS:

• To change the massage head, gently twist it and pull it out.

MAINTENANCE AND CARE:

- 1. Use a slightly damp towel to clean device and a dry with a soft cloth.
- 2. For storage or travel, remove battery and store in supplied carton box.

DOS & DON'TS:

- 1. Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation ports free from dust and debris.
- 2. Do not remove screws or attempt to disassemble.
- 3. Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.
- 4. Unplug the charger from device battery after charging or prior to use.
- 5. For adult use only. Do not use if injured. Consult your doctor before using this product.

Terms and Conditions:

- 1. This warranty is void if:
- a. The completely filled warranty card is not presented at the time of servicing the product.
- b. The product is not operated according to the instructions given under the user manual.
- c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
- d. Product has been damaged due to installation. epairs, alterations or modifications by unauthorized service organizations or persons.
- e. Product label specifying the model number, serial number and production code has been removed and altered.
- f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains. scratches, dents on the casing or paintwork of the product.
- g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
- h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
- 2. Repairs or replacements will be carried out by authorized service provider.
- 3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
- 4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be borne by the consumer.

- 5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warrantied, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the lass.
- 6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
- 7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
- 8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Consumer Warranty Card

Dear Customer.

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high au pro the

providing not only sen	as well as easy installation and use. At Lifelong, we believe in vice, but adding value to your purchase. The warranty has ned especially for you with your interest at heart.
them date of purchase Name of the product:	
	e receipt is necessary for warranty verification. omercare@lifelongindia.com
Address: Home Number: Office Number:	
Serial No.: ————————————————————————————————————	
	elongindiaonline.com and complete the online warranty form with you

personal and product details within 14 days.