

INTRODUCTION:

Dear Customer,

Congratulations! We welcome you to the Lifelong family, as you welcome our Lifelong Rechargeable Electric Toothbrush to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.


We are sure you can't wait to start using your new Lifelong Rechargeable Toothbrush.

We know we can't!

Looking forward to building our bond for life, cheers!

Thank you, Team Lifelong

1

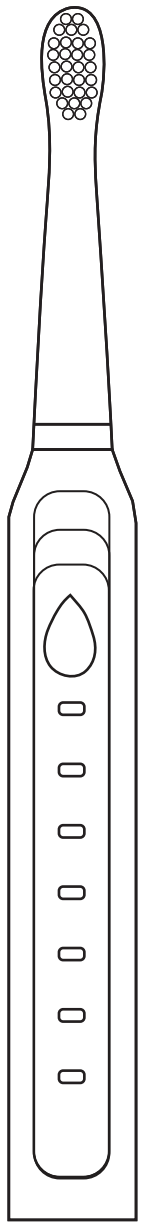

LLDC63

2

TABLE OF CONTENTS:


1. About the product
2. Marking of Product
3. Precautions Before Use
4. Operating Instructions
5. Product Features
6. Technical Features
7. Terms and conditions
8. Warranty

3

ABOUT THE PRODUCT:

Lifelong's Rechargeable Electric Toothbrush offers Sonic technology combined with strong brush action that gently removes plaque. Our unique technology brings you powerful, yet gentle clean, Powerful bristle vibrations give you a refreshing start every morning.

MARKING OF PRODUCT:

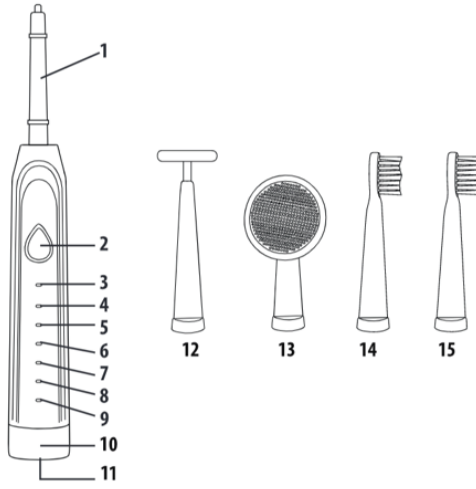

1. Mast for the interchangeable head
2. On/Off switch for work mode selection
3. Whitening mode
4. Massage mode
5. Cleaning mode
6. Sensitive teeth and gums mode
7. Polished mode

4

8. Gum care mode
9. LED charging indicator
10. Waterproof cover
11. USB charging port
12. Facelift bar head
13. Silicone brush head
14. Dentalcare brush head
15. Maxclean brush head

PRECAUTIONS BEFORE USE:

- Please do not use it outside the specified use.
- When connected with power, please check whether the power voltage is suited for your local voltage or not.
- Please use the original charger to avoid dangerously.
- Please avoid abrasives with whitening and removing tooth stains to use together.
- Please do not high strength use for teeth and gums.
- Children and can not operate their own or the mouthfeel weaker, please do not use.
- If the teeth are being treated, or If the Implanted tooth is not strong, do not use it.
- Do not drop It, do not hit It violently.
- Do not soak in water for long periods
- This product is only for personal use and is not available to dental clinical or organizations for multiple patients.

5

- Children and persons with physical and mental disabilities are required to use under the supervision of the guardian in the case of periodontal disease and oral symptoms, please use under the guidance of a doctor.
- If you feel uncomfortable or hurt in your mouth, please stop using it and consult a dentist. When the electric toothbrush fails, please consult your local government on how to dispose of this tool to maximize the protection of the environment.

OPERATING INSTRUCTIONS:

- Wet the toothbrush head and squeeze appropriate toothpaste on it.
- Place the bristles on the teeth of the 45-degree angle aligned with the gum line.
- Press the ON/OFF button, brush starts working/stops working.
- The product has 6 modes, choose any mode that suits you by pressing the power button successively to change the modes.

6

- You just gently move the bristles to place where need to clean, the pressure should be within 150g (Sonic toothbrush is to rely on high-frequency vibration of the bristles to clean the mouth, so when brushing teeth, the using force teeth is much smaller than ordinary, about 1/8 Orso can)
- The dentist advice that the oral cavity is divided into 4 areas (as shown above), slowly brushing the toothbrush during the brushing process to ensure that each area has about 30 seconds of cleaning time, the machine has a 30-second area reminder function, every 30 Seconds will be suspended for 0.5 seconds to remind you to replace the area to brush your teeth.

FACIAL CLEANSING:

1. Lather your face with soap and use the silicone brush attachment to thoroughly cleanse your skin.
2. Turn the power button to turn the device on and choose the mode you find comfortable for your face.

FACE LIFT:

1. Apply some facial cream onto your face and using the facelift bar attachment on your sonic brush body, gently massage your face going upwards.
2. Turn the power button on to use the vibration modes as well.

7

PRODUCT FEATURES:

- Ergonomic design: waterproof non-slip, feel comfortable, safe, and reliable
- USB fast-charge design, Unplug the waterproof plug at the bottom of the toothbrush. Plug the DC plug of the charging cable into the charging port at the bottom of the tooth brush before inserting it, and then insert the USB plug into the computer or the USB plug on the 5V.
- 2 brush heads & 2 Massage heads: standard with two soft brush heads and two facial massage heads, which can be replaced anytime, anywhere.

TECHNICAL FEATURES:

Name	Rechargeable Electric Toothbrush
Model Number	LLDC63
Rated Voltage	5V
Rated Power	1.4W

8

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

9

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

10

Consumer Warranty Card

Dear Customer,
Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from the date of purchase.

Name of the product: _____
Model: _____
Warranty coverage: 1 year
Please note:
Purchase receipt is necessary for warranty verification.
Customer Care: customercare@lifelongindia.com

Customer Details

Name: _____
Address: _____
Home Number: _____
Office Number: _____
E-mail Address: _____

Product Details

Model No.: _____
Serial No.: _____
Purchase Date: _____
Invoice Number: _____
Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.

11