

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

For queries and complaints:  
please contact: [customercare@lifelongonline.com](mailto:customercare@lifelongonline.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)

lifelong

# GUN MASSAGER

User Manual



**Built for Relaxation. Built for India.**

## INTRODUCTION:

Dear Customer,

Congratulations! We welcome you to the Lifelong family, as you welcome the Lifelong Gun Massager to be a part of your home.

We started Lifelong Online with a simple dream of being a part of every Indian household, which is why we have produced a varied assortment of innovative products across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle, and Smart Home - to better assist you.


We hope you enjoy using our Gun Massager as much as we enjoyed designing it for you.

Delivering promises that last lifelong.

Thank you for choosing us,

Lifelong

Built For India

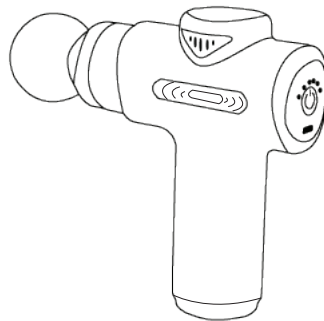

**LLGM109**

## TABLE OF CONTENTS:

1. About the product
2. Package Contents
3. Charging Instructions
4. Operation Instructions
5. Maintenance & Care
6. Dos and Don'ts
7. Terms & Conditions
8. Warranty Card

## ABOUT THE PRODUCT:

Lifelong Gun massager is a practical handheld massager that can relieve the body's soreness through deep and powerful muscle massage. It restores the body tissue, promotes muscle repair, and enhances the blood flow with five signature speed adjustments. It includes four ergonomic massage heads that cater to every muscle group for a comfortable and effective massage experience.

## PACKAGE CONTENTS:

- Massager
- Charger
- Seven Massage Heads: Flat Head, Spherical, U-shaped, Bullet head, Finger, Cylinder head, Shovel Head
- Manual
- Case

## MASSAGE HEAD:

1. Spherical Head: Arm, waist, hip, thigh, calf, and other large muscle groups
2. Flat Head: Neck, spine, and Achilles tendon
3. U-shaped Head: Deep tissue such as neck, spine, and Archilles tendon.
4. Bullet head: Relaxes deep tissues, such as meridians, palms, soles, etc.

- 5. Cylindrical: For precise deep massage around joints
- 6. Shovel Head: For relaxation and Ab toning
- 7. Finger: For acupressure massages

## CHARGING INSTRUCTIONS:

- Power button (on/off).
- The user can change the speed level when press the power button.
- Easy to replace massage heads by taking the old one out lightly.
- Recharging it with 5V output power adapter.
- Please recharge it about 3 hours before using at first time through the power adapter (5V output).
- For button model only: the 6 led lights will flash alternately when charging, when all the LED lights are on and solid/none are flashing, the massager is fully charged. when massaging, the LED light shows the current speed level.

## OPERATING INSTRUCTIONS:

- Power on the device by pressing the power button.
- With power switch in ON position, press the power button to adjust the speed level suitable for you.

- Press the power button again for switching off the device.
- Massage desired body part at the desired speed applying pressure as needed.

## MAINTENANCE AND CARE:

- To keep this device in its pristine condition, wipe its surface with a slightly damp towel and dry it with a soft cloth.
- Keep the device away from liquids and heat.

## DOS AND DON'TS:

1. This product should only be used on dry, clean skin, not through clothing. Please press and move gently on the skin. The duration of each section is about 60 seconds.
2. Do not use the device on the head or hard or bony parts of the body, parts with dense lymph, nerves and blood vessels, internal organs.
3. Please choose the right massage head for the body part you want to massage.
4. Bruising may occur with moderate or even pressure, always check the area being massaged.
5. Please adjust the equipment's head according to your body parts and pain. If you feel uncomfortable, please stop using it immediately.
6. Fingers, hair, or other body parts should be kept away from the shaft and head of the vibrator as they may be pinched.

7. Do not place anything in the ventilation openings of the device.
8. Do not immerse it in water or allow water to enter the vents of the device.
9. Do not drop or misuse the device, only charge it with a 5V DC Power adapter.
10. Please check the device and battery carefully before each use.
11. Do not alter or modify the device in any way.
12. Do not operate or charge unattended.

#### Notes:

- The interface of this product needs to be stretched and massaged. To prevent wear and tear, lubricating oil is applied. Be careful not to contaminate the clothes. If you accidentally touch the skin, you can wipe or wash it. It will disappear after a few uses.
- This product cannot be used while charging. The button model needs to long press the key to turn on and off.

#### Terms and Conditions:

1. This warranty is void if:
  - a. The completely filled warranty card is not presented at the time of servicing the product.
  - b. The product is not operated according to the instructions given under the user manual.
  - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
  - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
  - e. Product label specifying the model number, serial number and production code has been removed and altered.
  - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
  - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
  - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

## Consumer Warranty Card

Dear Customer,  
Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongonline.com](mailto:customercare@lifelongonline.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log on to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.