

Manufactured & Marketed by :  
**Lifelong Online Retail Private Limited**

For queries and complaints:  
please contact: [customercare@lifelongonline.com](mailto:customercare@lifelongonline.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)

# lifelong

## FOOT & CALF MASSAGER

User Manual



**Built for Relaxation. Built for India.**

## INTRODUCTION:


Dear Customer,

Congratulations!

Welcome to the world of Lifelong! Just like you're bringing our Lifelong Foot & Calf Massager into your life, we've made it with care to fit your daily needs. Here at Lifelong, we meticulously design our products to cater to your everyday requirements. This is the reason behind our diverse product lineup aimed at enhancing your comfort. We subsequently offer these products at highly competitive prices, with the aim of providing you and your family a soothing experience throughout the day.

We know you're excited to use your new Foot & Calf Massager. We're excited too, and we look forward to being friends for a long time. Cheers to the exciting journey ahead!

Thank you for choosing us,  
Lifelong  
Built For India

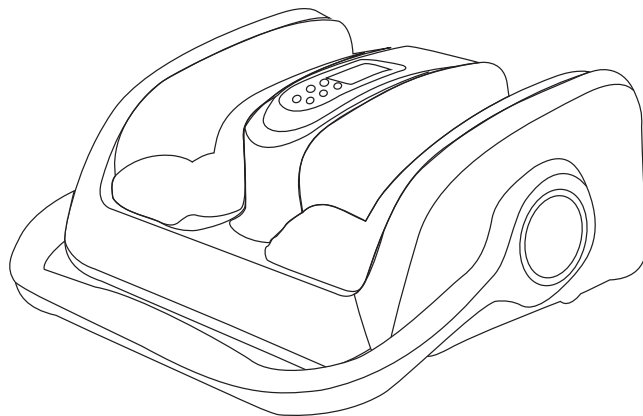

**LLM873**


## TABLE OF CONTENTS:

1. About the Product
2. Product Structure
3. Operating Instructions
4. Maintenance & Storage Method
5. Safety Precautions
6. Terms and Conditions
7. Warranty Card

## ABOUT THE PRODUCT:

Featuring a variety of massage techniques, including kneading, rolling, vibration, air compression, and heating, the Lifelong Foot & Calf Massager comes with a **360° rotatable handlebar** and offers an incredible experience. With the help of the handlebar, you can point the massager in the direction required and the part of the leg you want to relax. It also relieves tension, fatigue, improves circulation and sleep quality.

## PRODUCT STRUCTURE:

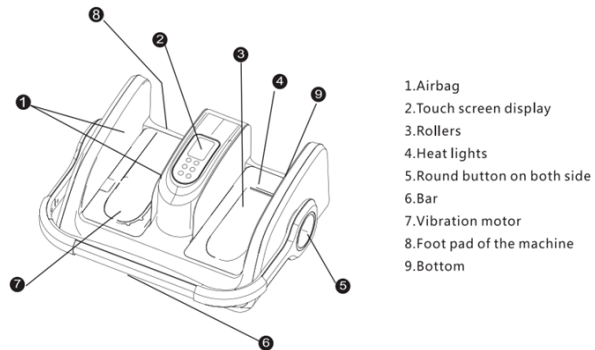

## OPERATING INSTRUCTIONS:



- Turn ON/OFF the machine – Start the machine at “15” mins, “F3” air pressure level, and “P1” auto program.



- There are five levels to choose. Choose the air pressure strength from F1 to F5.



- The Auto program of the machine gives an option to choose the cycle between: “P1”, “P2”, and “P3”.



- The Time function of the machine can be used to set the duration of the massage from 5 minutes to 30 minutes at an interval of every 5 minutes.



- To operate the Heat function, turn it on or off as per your choice.



- To massage the foot, ankle, or upper part of the calf, turn on/off the vibration function.

- To adjust the side handlebar, press the round button on both ends of the bar together and raise the machine. Release the button when the machine reaches the right angles.

## MAINTENANCE AND STORAGE METHOD:

- Ensure that you remove the plug from the electrical outlet before cleaning the massager.
- Do not use benzene, thinner or other abrasive solvents on your massager.
- Use a dry cloth while cleaning the massager, remember not to use a damp cloth to clean it.
- If needed, detach the fabric covers to clean and ensure good hygiene.
- When the product is not used for a long time, it is advisable to coil the power cord of the massager and store in a dust and moisture free location.
- Do not store the massager near intense heat or an open flame. Keep it away from direct sunshine.

## SAFETY PRECAUTIONS:

### Usage Instructions:

1. Do not use the massager in an excessively humid environment or around dust.
2. Do not use the massager in a hot room over 40°C.
3. Use the massager indoor only.
4. For those seeking a medical treatment or undergoing an ailment should consult their doctor before using the massager.
5. Other conditions in which people should avoid using the product are cancer, osteoporosis, high fever, impaired consciousness, diabetes, gout, malignant tumour, tuberculosis, severe cardiac dysfunction, bleeding disorders or tendencies, pregnancy, and other allergies.

## Warning:

1. Keep the massager away from children or people who cannot operate the product themselves.
2. Do not use the massager if the fabric cover is torn or missing.
3. Do not cover the massager while using it, as it can result in the heating up of the motor and subsequent breakdown or fire.
4. Do not stand or sit on the massager.
5. Use the parts given only by the manufacturer.
6. Do not insert or keep any object in the opening of the product.
7. Do not use the massager while sleeping or vice versa.
- 8. Do not use this product as a medical device.**
9. Avoid using the massager if you are intoxicated or unwell.
10. Use the massager an hour before or after eating.
11. Remember to turn off the massager after 30 minutes maximum
12. Refrain from using the massager if water is spilled on it.
13. If you feel any abnormal discomfort while using the massager, consult your doctor.
14. Do not use the massager if you experience malfunctioning in it.
15. Avoid using the product during a power failure or lightening.

## Terms and Conditions:

1. This warranty is void if:
  - a. The completely filled warranty card is not presented at the time of servicing the product.
  - b. The product is not operated according to the instructions given under the user manual.
  - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
  - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
  - e. Product label specifying the model number, serial number and production code has been removed and altered.
  - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
  - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
  - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

## Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongonline.com](mailto:customercare@lifelongonline.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log on to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.