

Manufactured & Marketed by :
Lifelong Online Retail Private Limited

For queries and complaints:
please contact: customercare@lifelongonline.com
www.lifelongindiaonline.com

lifelong

AIR PRESSURE MASSAGER

User Manual



Built for Relaxation. Built for India.

INTRODUCTION:


Dear Customer,

Congratulations!

Welcome to the world of Lifelong! Just like you're bringing our Lifelong Air Pressure Calf Massager into your life, we've made it with care to fit your daily needs. Here at Lifelong, we meticulously design our products to cater to your everyday requirements. This is the reason behind our diverse product lineup aimed at enhancing your comfort. We subsequently offer these products at highly competitive prices, with the aim of providing you and your family a soothing experience throughout the day.

We know you're excited to use your new Lifelong Air Pressure Calf Massager. We're excited too, and we look forward to being friends for a long time. Cheers to the exciting journey ahead!

Thank you for choosing us,
Lifelong
Built For India

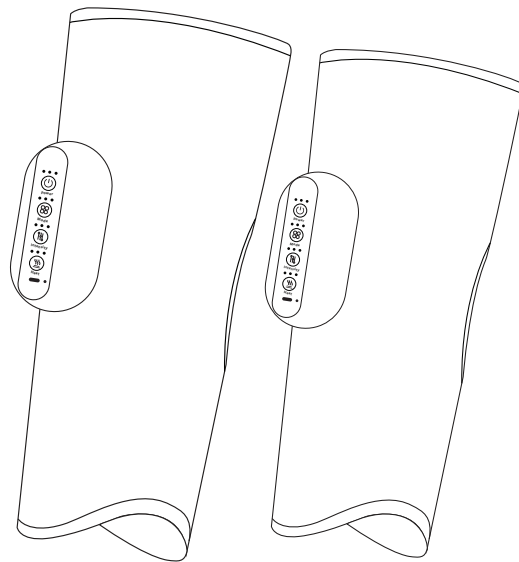

LLM882 / LLM864

TABLE OF CONTENTS:

1. About the Product
2. Package Contents & Name of Components
3. Usage Instructions
4. Maintenance & Storage Method
5. Safety Precautions
6. Troubleshooting
7. Terms and Conditions
8. Warranty Card

ABOUT THE PRODUCT:

The rechargeable Lifelong Air Pressure Calf Massager is here to provide the ultimate relief. This innovative device uses gentle air pressure to massage your calf muscles, promoting improved circulation and easing away tension.

PACKAGE CONTENTS & NAMES OF COMPONENTS:


Package Contents:

Once you unpack, check the package contents before use. The package should have:

- 2 Massagers
- USB Power Line
- Instruction Manual


Component Names:

Front of Massager:

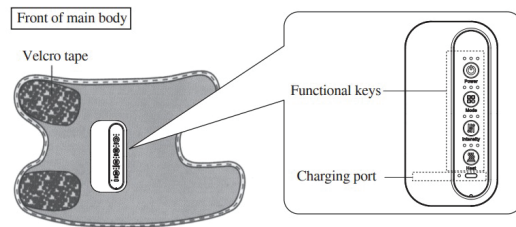

Back of Massager:


Back of main body

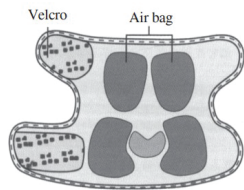

USAGE INSTRUCTIONS:

Preparation before use:

Connect the massager with the charging line.

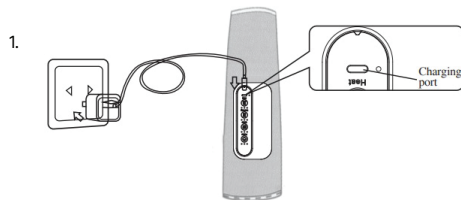

1. Insert the Type C end of the USB cable into the charging port of the massager.

- Insert the Type A end of the USB cable into the USB adapter port.
- Ensure that the USB cable is properly inserted to the bottom.


2. Start Charging

- The charging begins when the indicator light flashes. Once the charging is complete, the indicator light is switched on.
- The indicator light flashes when the light is low.
- The battery takes 1-2 hours to get fully charged, and it can be used till 90 minutes once it is charged.

Additional Note:

- In case you are charging the device at a low ambient temperature, automatic protection may get initiated due to which charging may be paused.
- When the product is used for the first time or not used for a long time, recharge the battery fully.
- Do not use the product while it is getting charged.


DESCRIPTION OF BUTTONS:

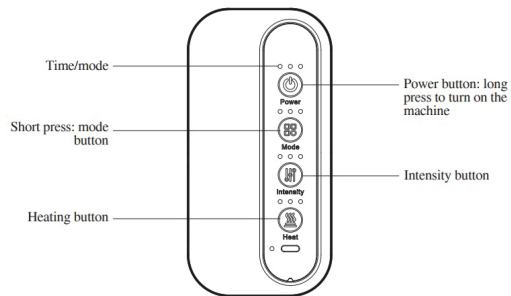

Start Using the Product:


1. Wear the Massager

Adjust the position of the main body and fix it with the rear Velcro tape.

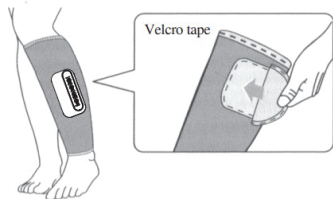


2. Long Press the Power Button

Long press the power button for 3 seconds, and the machine gets turned on/off with a beep. The machine starts with mode (1), intensity (1), heating (1), and the indicator light gets switched on. Short press the power button to adjust the time. There are three adjustable time modes.

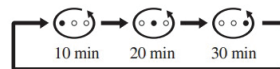


3. Adjustment


Mode: Short press the mode key to switch to different massage modes from mode 1, mode 2 and mode 3, displayed on the pilot lamp.

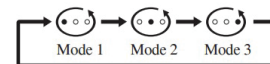

Intensity: Short press the intensity button to switch to the desired intensity. (Intensity 1, Intensity 2, and Intensity 3), displayed on the pilot lamp.

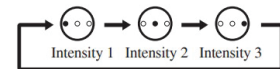

Temperature: Short press the temperature button to switch to different temperatures. (1 for low, 2 for medium, 3 for high, and turn off)

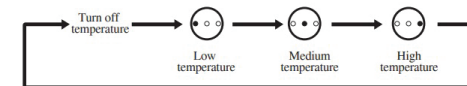

MAINTENANCE & STORAGE METHOD:


Maintenance Method:

1. Controller Care

Handle the screen of the machine with care, and use a soft, clean cloth to wipe it.

2. Cleaning the Fabric

The fabric cannot be washed. Use alcohol to kill bacteria or remove sweat or use a clean wet towel to wipe it.

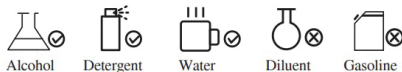

Warning

- Do not wash the main body, USB line or remote controller. Otherwise, it may cause electric shock, short-circuit or fire.

- Keep the product in a dry place with low humidity after cleaning it from stains and dust.
- If the product is not used for a long time, cover it with a cloth and keep it in a dust-free place.
- Keep the product out of the reach of babies or children.
- Arrange the power cord neatly before putting it away.

SAFETY PRECAUTIONS:

The consumers should use the product prudently or under the guidance of a physician

- A user with leg nerve disorder, skin perception disorder, allergic constitution, skin disease or odorous dermatitis should not use the product.
- A user who suffers or had suffered from leg disease (i.e., tendon injury, arthritic deformity, fracture, and varicose veins) should avoid using this product.
- A user who is sensitive to heat should avoid the product.
- A person whose pacemakers are susceptible to electromagnetic failures should not use the product.
- A person with heart disease should not use the product.
- A woman in an unstable pregnancy condition or one about to give birth should not use the product.
- Do not operate the massager under a pillow or a blanket to avoid fire, electric shock, or injury due to overheating.
- Do not use the product in a high temperature environment or humidity, to avoid electric shock.
- **Do not use the product for medical purposes.**
- Do not disassemble the internal components of the product without authorization.
- Do not use the product while sleeping.
- Do not use the massager if you are intoxicated or feel discomfort.

TROUBLESHOOTING:

In case of any issue in using the product, refer to the troubleshooting steps given below.

If the product still fails to work, contact the nearest service center.

Phenomena	Possible Reasons	Treatment
Abnormal sound during operation	Noise during normal operation or noise due to friction between internal components and cloth cover	No treatment is needed
No operation or no perceptible intensity	Inappropriate tightness The product has not started Low battery	Re-attach the Velcro Tape Start the product Recharge the battery
Shutdown during use	Times up by default	Re-start after 15 minutes

Terms and Conditions:

1. This warranty is void if:
 - a. The completely filled warranty card is not presented at the time of servicing the product.
 - b. The product is not operated according to the instructions given under the user manual.
 - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
 - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
 - e. Product label specifying the model number, serial number and production code has been removed and altered.
 - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
 - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
 - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

Warranty Service

All Lifelong consumer products are covered against manufacturing defects from their date of purchase.

Name of the product: _____

Model: _____

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: customercare@lifelongonline.com

Customer Details

Name: _____

Address: _____

Home Number: _____

Office Number: _____

E-mail Address: _____

Product Details

Model No.: _____

Serial No.: _____

Purchase Date: _____

Invoice Number: _____

Online Site: _____

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.