

RECHARGEABLE HEAD, SCALP AND FULL BODY PAIN RELIEF MASSAGER

User Manual



Built for Relaxation. Built for India.

Manufactured for & Marketed by :

Lifelong Online Retail Private Limited

For queries and complaints: please contact: customercare@lifelongindia.com www.lifelongindiaonline.com

INTRODUCTION:

Dear Customer, Congratulations!

We welcome you to the Lifelong family, as you welcome our Lifelong

Head, Scalp and Full Body Massager to be a part of your home.

At Lifelong, our products are carefully built for the needs of India, which is why we

crafted a diverse range of innovative products, across a spectrum of categories -

Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you.

We then deliver them at the best prices, enriching your everyday living.


We are sure you can't wait to start using your new Lifelong Head, Scalp and

Full Body Massager.

We know we can't!

Looking forward to building our bond for life, cheers!

Thank you, Team Lifelong

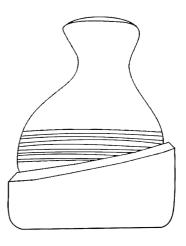

LLM765

1


TABLE OF CONTENTS:

- 1. About the product
- 2. Product Markings
- 3. Safety Precautions
- 4. Using the Product
- 5. Charging
- 6. Device Cleaning & Maintenance
- 7. Terms & Conditions
- 8. Warranty Card

ABOUT THE PRODUCT:

The Lifelong Head, Scalp and Full Body Massager reduces pain and promotes blood circulation. The massager increases blood circulation to the scalp, promoting relaxation, releasing stress and boosting hair growth. The nodes are waterproof and removable. Being cordless and rechargeable, this massager can be carried anywhere so that you enjoy the perfect massage anywhere without an annoying wire creating any hassles. For hair loss in the head, less hair, poor hair quality, the massager has a significant effect and can significantly improve fatigue.

PRODUCT MARKINGS:

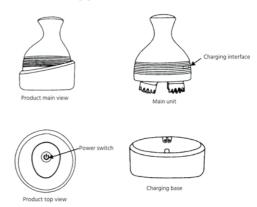

SAFETY PRECAUTIONS:

- Do not put any heavy weight on the massager.
- Clean the device with a slightly damp cloth.
- · Do not wash any part of the massager.
- Do not disassemble the instrument by yourself.
- Patients with severe stroke, hypertension should use the massager with caution.
 Please consult a professional physician. If there is any discomfort during use, stop using it immediately.
- Please avoid using the massager after drinking or when you are not feeling well.
- Please use the power adapter configured with this product and pay attention to electricity safety.
- The massager is designed as a home electronic health care product and is not suitable for diagnostic and therapeutic purposes.
- Do not place the charging base in a humid environment.
- Do not charge in excessively humid conditions.

Note:

- · This product is for household use only
- · Specifications are subject to change without prior notice

USING THE PRODUCT:

1. Press the power button to switch on the massager. The power button is also used to manage the intensity and speed of the massager.

- 2. There are four speed modes:
- Mode 1: The default speed is slow.
- Mode 2: The second speed is fast.
- Mode 3: The third speed is anti-clockwise slow.
- Mode 4: The fourth is anti-clockwise fast.
- 3. Press and hold for 3 seconds to switch off.
- 4. The massager comes with a power adapter.
- 5. If the machine battery is too low, the battery indicator light turns red, the machine should be charged to use again.

CHARGING:

- During the charging process, the battery indicator light is always red.
 When the battery is fully charged, the indicator light turns vellow.
- 2. Do not use this product during charging (the charging is not supported in the power-on state. The product starts the self-protection program and automatically
- stops working.
- 3. The charging time is 2 hours.

DEVICE CLEANING & MAINTENANCE:

- Moisten a soft cloth with medical alcohol or clean water to wipe the machine's body to remove the dirt.
- Please keep this machine in a safe place after use. No heavy weight is allowed on the machine.
- Please ask your professional technician for repair if any fault occurs.

Terms and Conditions:

1.This warranty is void if:

- a. The completely filled warranty card is not presented at the time of servicing the product.
- b. The product is not operated according to the instructions given under the user manual.
- C. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
- d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
- e. Product label specifying the model number, serial number and production code has been removed and altered.
- Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
- g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
- Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
- Repairs or replacements will be carried out by authorized service provider.
- During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelona and/or its authorised servicer.
- All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be

- 5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warrantied, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the lass.
- 6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
- 7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profts, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
- No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

8

Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

date of purchase. Name of the prod	mer products are covered against manufacturing defects from them
	nase receipt is necessary for warranty verification. Istomercare@lifelongindia.com
Address: Home Number: Office Number:	ls
Product Details Model No.: —— Serial No.: —— Purchase Date: Invoice Number:	

Please log on to www.lifelongindiaonline.com and complete the online warranty form with your personal and product details within 14 days.