

# lifelong

## MINI MASSAGER

User Manual



Built for Relaxation. Built for India.

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

For queries and complaints:  
please contact: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)

## INTRODUCTION:

Dear Customer,

Congratulations! We welcome you to the Lifelong family, as you welcome our Lifelong Mini Head and Body Massager to be a part of your home.


At Lifelong, our products are carefully built for the needs of India, which is why we have crafted a diverse range of innovative products, across a spectrum of categories - Home, Kitchen, Grooming, Fitness, Lifestyle and Smart Home, to better help you. We then deliver them at the best prices, enriching your everyday living.

We are sure you can't wait to start using your new Lifelong Mini Head and Body Massager.

We know we can't!

Looking forward to building our bond for life, cheers!

Thank you, Team Lifelong

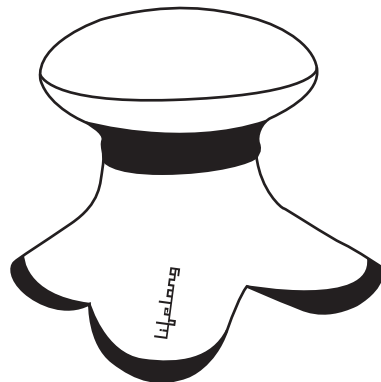

**LLM774**

## TABLE OF CONTENTS:

1. About the product
2. Additional Features
3. Battery Installation
4. How to use
5. Device Cleaning & Maintenance
6. Safety for usage/ Dos and Don'ts
7. Product specifications
8. Terms and Conditions
9. Warranty Card

## ABOUT THE PRODUCT:

The Mini Head and Body Massager is a rejuvenating device that gives you great value for money. The massager kneads and relaxes your muscles, letting the days worries melt away. With more than 100 innovative products available , no wonder millions of customers continue to choose Lifelong! It has strongly been trusted by loyal users for over five years. Made for you, and made especially for your relaxation.

## ADDITIONAL FEATURES:

- Enjoy effortless relaxing massage anywhere, any time with this vibrating mini massager.
- Its ergonomic design and its 3 robotic massage contacts allow a proper massage all along your body. Also, its vibrating movement benefits in blood circulation and helps in muscle relaxation.
- When combined with a balanced diet and exercise, the massage helps to decrease fat deposits.

## BATTERY INSTALLATION:



1. Main Unit



2. Batteries needed  
(AAAx3)



3. User Manual

1. Open the cover of the main machine according to the indicated direction and put 3 AAA batteries inside. (Battery not included)

2. Close the battery cover

- Please confirm the anode and cathode of batteries.

## HOW TO USE:

1. Place the 3 massage fingers of the device on the part you wish to massage, and press the top of the main machine with proper force for it to start working.

2. Detach the device from your body and release the button on the massager to stop the massage.

## DEVICE CLEANING & MAINTENANCE:

- Moisten a soft cloth with medical alcohol or clean water to wipe the machine's body to remove the dirt.
- Please keep this machine in a safe place after use. No heavy weight is allowed on the machine.
- Please ask your professional technician for repair if any fault occurs.

## DOS & DON'TS:

- Do not put any heavy weight on the massager.
- Clean the device with a slightly damp cloth.
- Do not wash any part of the massager.
- Do not disassemble the instrument by yourself.

## PRODUCT SPECIFICATIONS:

- Dimension: 105.6mm x 6.4mm x 89.8mm
- Battery: AAA dry battery (Battery not included)
- Net Weight: 204.5 g
- Materials: ABS
- Power consumption: 1.3 VA
- Ingression Protection Rating: IP67

## Terms and Conditions:

1. This warranty is void if:
  - a. The completely filled warranty card is not presented at the time of servicing the product.
  - b. The product is not operated according to the instructions given under the user manual.
  - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
  - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
  - e. Product label specifying the model number, serial number and production code has been removed and altered.
  - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
  - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
  - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as will as any other expenses and incidentals will be borne by the consumer.

5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not reply on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

## Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log on to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.