

## Terms and Conditions:

1. This warranty is void if:
  - a. The completely filled warranty card is not presented at the time of servicing the product.
  - b. The product is not operated according to the instructions given under the user manual.
  - c. Damages are caused by lightning, abnormal voltage, water or other liquid intrusion, fire, flood, accident, negligence or improper handling.
  - d. Product has been damaged due to installation, repairs, alterations or modifications by unauthorized service organizations or persons.
  - e. Product label specifying the model number, serial number and production code has been removed and altered.
  - f. Defects or parts requiring replacement due to ordinary wear and tear, corrosion, rust or stains, scratches, dents on the casing or paintwork of the product.
  - g. Claims damaged and/or missing parts (accessories) after 7 days from the original date of goods received.
  - h. Defects or faults in product which have been used for commercial/industrial purposes or which have been rented/leased or have been otherwise subject to non-household/non-domestic use.
2. Repairs or replacements will be carried out by authorized service provider.
3. During the limited warranty period, Lifelong or its authorized service provider will repair without charge the defective unit inclusive of labor and parts and restore the unit to its optimum working condition. All defective parts used for the warranty repair should be surrendered to Lifelong and/or its authorised servicer.
4. All expenses incurred in collecting the unit (s) or part(s) thereof from authorised service provider as well as any other expenses and incidentals will be borne by the consumer.
5. Lifelong obligations are limited to the repair and replacement of the defective product. Except as set forth above, there are no other express or implied warranty and all warranted, conditions or other terms implied by statute or common law (including any warranty of satisfactory quality, merchantability or fitness for a particular purpose) are excluded to the fullest extent permitted by the law.
6. Lifelong total liability for damages relating to or arising out of the purchase or use of the product regardless of the type or cause of such damage of the form of characterization of the claim asserted (e.g. contract or tort) shall not exceed the original purchase price paid for the product.
7. However in no event shall Lifelong and Lifelong authorized distributors be liable for any punitive, special incidental, indirect or consequential losses or damages whatsoever (including without limitation, damages for the loss of revenue, business, profits, goodwill, or contracts, business interruptions, loss of business information or any other pecuniary loss.) Whether or not Lifelong has been advised of the possibility of such damages. These limitations shall apply not with standing the failure of the essential purpose of any limited warranty. This limited warranty does not affect consumer's statutory rights under the law.
8. No carrier, dealer or employee is thereof authorized to make modifications to this warranty and you should not rely on any such representation. Lifelong reserves the right to amend the terms and conditions if necessary.

Manufactured for & Marketed by :  
**Lifelong Online Retail Private Limited**

Made in China

For queries and complaints:  
please contact: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)  
[www.lifelongindiaonline.com](http://www.lifelongindiaonline.com)

# lifelong

## Full Body Massage Chair with Recliner LLM549



POWERFUL 3D BACK  
& LEG MASSAGE



IN-BUILT BLUETOOTH  
SPEAKERS



TARGET SPECIFIC  
ZONES



ZERO GRAVITY  
RECLINE



HEATING AND  
VIBRATION



## Consumer Warranty Card

Dear Customer,

Thank you for choosing a Lifelong consumer product. All Lifelong Consumer products are designed and manufactured to the highest standards to deliver high quality performance, as well as easy installation and use. At Lifelong, we believe in providing not only service, but adding value to your purchase. The warranty has therefore been designed especially for you with your interest at heart.

### Warranty Service

All Lifelong consumer products are covered against manufacturing defects from them date of purchase.

Name of the product: \_\_\_\_\_

Model: \_\_\_\_\_

Warranty coverage: 1 year

Please note: Purchase receipt is necessary for warranty verification.

Customer Care: [customercare@lifelongindia.com](mailto:customercare@lifelongindia.com)

### Customer Details

Name: \_\_\_\_\_

Address: \_\_\_\_\_

Home Number: \_\_\_\_\_

Office Number: \_\_\_\_\_

E-mail Address: \_\_\_\_\_

### Product Details

Model No.: \_\_\_\_\_

Serial No.: \_\_\_\_\_

Purchase Date: \_\_\_\_\_

Invoice Number: \_\_\_\_\_

Online Site: \_\_\_\_\_

Please log onto to [www.lifelongindiaonline.com](http://www.lifelongindiaonline.com) and complete the online warranty form with your personal and product details within 14 days.

## PRODUCT MAINTENANCE

1. Please turn off the power switch after usage.
2. Roll up the appliance power cord if not in use for a long time, and store in a dry and clean environment.
3. Do not store in high temperature or near the fire, do not keep under direct sunlight for a long time.
4. Please use dry cloth to do the cleaning, do not use thinner, benzene or alcohol on equipment.
5. Please use it intermittently, you should not let this product be in continuous operation.

## COMMON PHENOMENON PROCESSING

1. A click sound during operation, the normal sound of working operation.
2. If instrument stops operation, please press key to operate again.

## SHUT OFF THE POWER AND STOP MESSAGE

During massage, if you press on-of key, the massager will shut off immediately and the massage chair will reset automatically .

When the massager completes its cycle of 20 minutes , all the massager functions also will shut off immediately, but the backrest of the chair and leg lifter and seat do not get automatic reset but are in the original shape and state.

Dear user,

Thank you for purchasing the Lifelong massage chair. Please read the details in this pamphlet carefully and pay special attention to the important safety instructions to ensure proper use. Please keep this brochure for your consultation.

## SAFETY AND MAINTENANCE INSTRUCTIONS

### Warning!

1. Do not use this machine with unsuitable standard voltage.
  2. Be sure to use the plug with matching socket.
  3. Please use the massage chair according to the operating manual.
  4. Please do not use the massage chair outdoor.
  5. The massage chair is not intended to be used for persons with body disabilities, senses and nerve defect (including children) unless under super vision and guidance.
- If you suffer from the below conditions, we advise you not to use the Massage chair:
1. Osteoporosis Patients.
  2. With heart problems, and those who wear electronic medical devices such as pacemakers.
  3. Pregnant or Menstruating Women.
  4. Physical trauma or body surface injured patients.

## SECURITY

1. Please check the voltage meets the specifications of the product.
2. Do not use wet hands to remove the power plug.
3. Never use a damp cloth to clean this machine's power switch, plugs, etc.
4. Always consult a doctor if you suffer from chronic pain.

## FUNCTIONS AND EFFICIENCY

1. There are 5 types of simulation massage techniques

i) **Kneading** -- the rollers move in a circular pattern.

ii) **Rolling** -- the rollers move up and down the backrest.

iii) **Tapping** -- the rollers alternate pushing in and out from the backrest, simulating the popular "karate chop" massage technique.


iv) **Gripping** -- the chair includes devices that grip your legs or arms in a snug hold before releasing.

v) **Shiatsu** -- the rollers press against certain points on the your back (sometimes called acupoints) to relieve tension.

2. Air pressure massage function for hips, calves and feet.

## USAGE OF THE LIFELONG MASSAGE CHAIR


- |   |  |
|---|--|
| 1. Eliminate fatigue                                | 2. Promote blood circulation           |
| 3. Relieve Muscle stress                            | 4. Reduce stress, tension, and anxiety |
| 5. Boost immune system                              | 6. Improve posture                     |
| 7. Reduces pressure on nerves by aligning the spine |  |

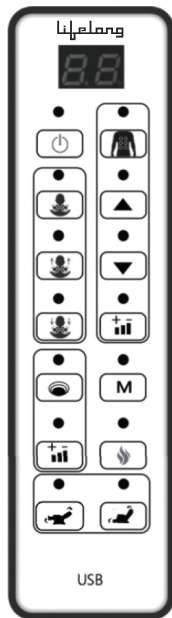

### ATTENTION

- In order to avoid the remote controller function failure, do not use the control panel with wet hand.
- In order to avoid damaging the remote controller, do not splash water or other liquids into the controller.



## 6. Back Intensity

Press  to set level of intensity for the back tapping function. Slow flashing of indicator light point to low level, fast flashing points to middle level, indicator light on points to high level.


## 7. Mode

The control panel has three modes and the indicator with three flashes which shows different tapping function.

## 8. Airbag

In the working state, press  to enjoy airbag function. The airbag massage is on both sides of the hips, middle of the hips, and the calf. Press  again to stop it.

## 9. Leg Intensity

Press  to set level of intensity for the leg function. Slow flashing of indicator light point to low level, fast flashing points to middle level, indicator light on points to high level.

## 10. Recliner

### • Down

The button takes the backrest down and the calf up.


### • UP

The button takes the backrest up and the calf down.


## 11. Bluetooth

Use the Bluetooth connect feature to enjoy music during the massage.

## NAMES OF PARTS AND FUNCTIONS

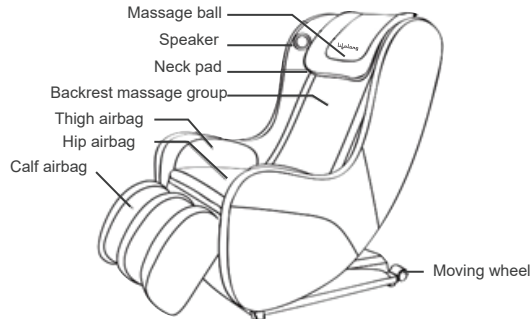


Full power switch diagram  
(Lower back of the massage chair)

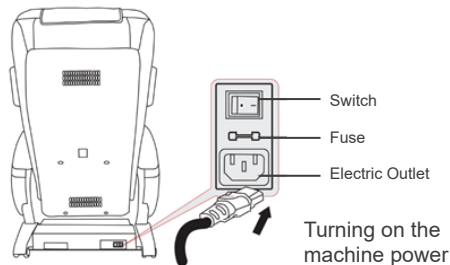

## INSTRUCTIONS OF CALF AND SEAT CUSHION

Turn the calf cushion down parallel to the seat

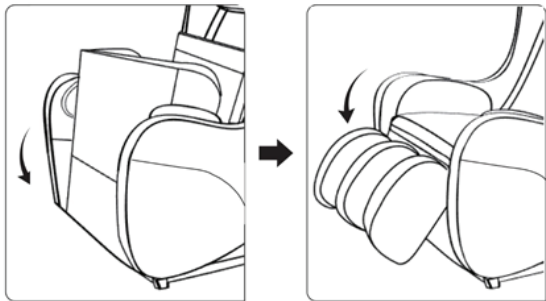


## MOVING CHAIR

Press the product back when moving the chair, in order to force moving roller. Move the chair forward or backward. Don't sit on the chair while moving.



## DISTRIBUTION DIAGRAM OF CONTROL PANEL

### 1. On/off

In standby mode, press  to start the massage chair. Chair will work for maximum 20 minutes (default time). Again press (icon) to turn it off.

### 2. Neck


In the working state, press (icon) to turn on or off the kneading function of the neck, and the corresponding indicator light will be turned on or off.

### 3. UP & DOWN





When the product is in standby or working state, long press the up button to make the neck massage ball move upward. Loosen the button, and the massage ball will stop moving upward.

Again long press the down button to make the neck massage ball move downward. Loosen the button, and the massage ball will stop moving downward.

### 4. Full-back button

In the working state, press  (indicator lights off) to enjoy kneading and tapping on the back. Press again (indicator lights on), for fixed point massage.

### 5. Upper back & lower back

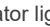In the working state, press  indicator lights on for kneading and tapping upper back. Press  again (indicator lights off) for fixed point upper back massage. Press , (indicator lights on) for kneading and tapping lower back massage. Press  again (indicator lights off) for fixed point lower back massage.